

Firmware documentation - Master
v. 6.1.5

Generated by Doxygen 1.8.13

Contents

1	Firmware	1
2	Data Structure Index	3
2.1	Data Structures	3
3	File Index	5
3.1	File List	5
4	Data Structure Documentation	7
4.1	st_calib Struct Reference	7
4.1.1	Detailed Description	7
4.1.2	Field Documentation	7
4.1.2.1	direction	7
4.1.2.2	enabled	7
4.1.2.3	repetitions	8
4.1.2.4	speed	8
4.2	st_data Struct Reference	8
4.2.1	Detailed Description	8
4.2.2	Field Documentation	8
4.2.2.1	buffer	8
4.2.2.2	ind	9
4.2.2.3	length	9
4.2.2.4	ready	9
4.3	st_dev Struct Reference	9
4.3.1	Detailed Description	9

4.3.2	Field Documentation	9
4.3.2.1	pwm_limit	10
4.3.2.2	tension	10
4.3.2.3	tension_conv_factor	10
4.3.2.4	tension_valid	10
4.4	st_meas Struct Reference	10
4.4.1	Field Documentation	10
4.4.1.1	acc	11
4.4.1.2	curr	11
4.4.1.3	emg	11
4.4.1.4	hand_meas	11
4.4.1.5	pos	11
4.4.1.6	rot	11
4.4.1.7	vel	11
4.5	st_mem Struct Reference	12
4.5.1	Field Documentation	12
4.5.1.1	activ	13
4.5.1.2	activate_pwm_rescaling	13
4.5.1.3	baud_rate	13
4.5.1.4	control_mode	13
4.5.1.5	curr_lookup	13
4.5.1.6	current_limit	13
4.5.1.7	double_encoder_on_off	13
4.5.1.8	emg_calibration_flag	13
4.5.1.9	emg_max_value	14
4.5.1.10	emg_speed	14
4.5.1.11	emg_threshold	14
4.5.1.12	flag	14
4.5.1.13	id	14
4.5.1.14	input_mode	14

4.5.1.15	k_d	14
4.5.1.16	k_d_c	14
4.5.1.17	k_d_c_dl	15
4.5.1.18	k_d_dl	15
4.5.1.19	k_i	15
4.5.1.20	k_i_c	15
4.5.1.21	k_i_c_dl	15
4.5.1.22	k_i_dl	15
4.5.1.23	k_p	15
4.5.1.24	k_p_c	15
4.5.1.25	k_p_c_dl	16
4.5.1.26	k_p_dl	16
4.5.1.27	m_mult	16
4.5.1.28	m_off	16
4.5.1.29	max_step_neg	16
4.5.1.30	max_step_pos	16
4.5.1.31	motor_handle_ratio	16
4.5.1.32	pos_lim_flag	16
4.5.1.33	pos_lim_inf	17
4.5.1.34	pos_lim_sup	17
4.5.1.35	res	17
4.5.1.36	rest_delay	17
4.5.1.37	rest_pos	17
4.5.1.38	rest_ratio	17
4.5.1.39	rest_vel	17
4.5.1.40	watchdog_period	18
4.6	st_ref Struct Reference	18
4.6.1	Detailed Description	18
4.6.2	Field Documentation	18
4.6.2.1	curr	18
4.6.2.2	onoff	18
4.6.2.3	pos	18
4.6.2.4	pwm	18

5 File Documentation	19
5.1 command_processing.c File Reference	19
5.1.1 Detailed Description	20
5.1.2 Function Documentation	21
5.1.2.1 cmd_activate()	21
5.1.2.2 cmd_get_activate()	21
5.1.2.3 cmd_get_currents()	21
5.1.2.4 cmd_get_emg()	21
5.1.2.5 cmd_get_inputs()	21
5.1.2.6 cmd_get_measurements()	21
5.1.2.7 cmd_get_measurements_from_SH()	21
5.1.2.8 cmd_ping()	22
5.1.2.9 cmd_set_baudrate()	22
5.1.2.10 cmd_set_inputs()	22
5.1.2.11 cmd_set_watchdog()	22
5.1.2.12 cmd_store_params()	22
5.1.2.13 commProcess()	22
5.1.2.14 commReadMeasFromAnother()	22
5.1.2.15 commWrite()	22
5.1.2.16 commWrite_old_id()	23
5.1.2.17 commWriteAnother()	23
5.1.2.18 get_param_list()	23
5.1.2.19 infoGet()	24
5.1.2.20 infoPrepare()	24
5.1.2.21 infoSend()	24
5.1.2.22 LCRChecksum()	24
5.1.2.23 memInit()	25
5.1.2.24 memRecall()	25
5.1.2.25 memRestore()	25
5.1.2.26 memStore()	25

5.1.2.27	sendAcknowledgment()	26
5.1.2.28	setZeros()	26
5.2	command_processing.h File Reference	26
5.2.1	Detailed Description	28
5.2.2	Function Documentation	29
5.2.2.1	cmd_activate()	29
5.2.2.2	cmd_get_activate()	29
5.2.2.3	cmd_get_currents()	29
5.2.2.4	cmd_get_emg()	29
5.2.2.5	cmd_get_inputs()	29
5.2.2.6	cmd_get_measurements()	29
5.2.2.7	cmd_get_measurements_from_SH()	29
5.2.2.8	cmd_ping()	30
5.2.2.9	cmd_set_baudrate()	30
5.2.2.10	cmd_set_inputs()	30
5.2.2.11	cmd_set_watchdog()	30
5.2.2.12	cmd_store_params()	30
5.2.2.13	commProcess()	30
5.2.2.14	commReadMeasFromAnother()	30
5.2.2.15	commWrite()	30
5.2.2.16	commWrite_old_id()	31
5.2.2.17	commWriteAnother()	31
5.2.2.18	get_param_list()	31
5.2.2.19	infoGet()	32
5.2.2.20	infoPrepare()	32
5.2.2.21	infoSend()	32
5.2.2.22	LCRChecksum()	32
5.2.2.23	memInit()	33
5.2.2.24	memRecall()	33
5.2.2.25	memRestore()	33

5.2.2.26	memStore()	33
5.2.2.27	sendAcknowledgment()	34
5.2.2.28	setZeros()	34
5.3	commands.h File Reference	34
5.3.1	Detailed Description	36
5.3.2	Enumeration Type Documentation	36
5.3.2.1	qbmove_command	37
5.3.2.2	qbmove_control_mode	38
5.3.2.3	qbmove_input_mode	38
5.3.2.4	qbmove_parameter	38
5.4	globals.c File Reference	39
5.4.1	Detailed Description	40
5.4.2	Variable Documentation	40
5.4.2.1	c_mem	41
5.4.2.2	dev_pwm_limit	41
5.4.2.3	dev_tension	41
5.4.2.4	dev_tension_f	41
5.4.2.5	forced_open	41
5.4.2.6	g_measOld	41
5.4.2.7	g_refOld	41
5.4.2.8	g_rx	41
5.4.2.9	interrupt_flag	42
5.4.2.10	master_mode	42
5.4.2.11	pow_tension	42
5.4.2.12	pwm_sign	42
5.4.2.13	receive_meas_from_hand	42
5.4.2.14	reset_last_value_flag	42
5.4.2.15	rest_enabled	42
5.4.2.16	rest_pos_curr_ref	42
5.4.2.17	tau_feedback	43

5.4.2.18	tension_valid	43
5.4.2.19	timer_value	43
5.4.2.20	timer_value0	43
5.4.2.21	watchdog_flag	43
5.5	globals.h File Reference	43
5.5.1	Detailed Description	46
5.5.2	Macro Definition Documentation	46
5.5.2.1	ANTI_WINDUP	46
5.5.2.2	CALIBRATION_DIV	47
5.5.2.3	CURR_INTEGRAL_SAT_LIMIT	47
5.5.2.4	CURRENT_HYSTERESIS	47
5.5.2.5	DEFAULT_CURRENT_LIMIT	47
5.5.2.6	DEFAULT_EEPROM_DISPLACEMENT	47
5.5.2.7	DIV_INIT_VALUE	47
5.5.2.8	EMG_SAMPLE_TO_DISCARD	47
5.5.2.9	LOOKUP_DIM	47
5.5.2.10	MAX_WATCHDOG_TIMER	48
5.5.2.11	NUM_OF_ANALOG_INPUTS	48
5.5.2.12	NUM_OF_EMGS	48
5.5.2.13	NUM_OF_MOTORS	48
5.5.2.14	NUM_OF_PARAMS	48
5.5.2.15	NUM_OF_SENSORS	48
5.5.2.16	POS_INTEGRAL_SAT_LIMIT	48
5.5.2.17	PWM_MAX_VALUE	48
5.5.2.18	RECEIVE	49
5.5.2.19	SAMPLES_FOR_EMG_MEAN	49
5.5.2.20	SAMPLES_FOR_MEAN	49
5.5.2.21	UNLOAD	49
5.5.2.22	WAIT_ID	49
5.5.2.23	WAIT_LENGTH	49

5.5.2.24	WAIT_START	49
5.5.3	Enumeration Type Documentation	49
5.5.3.1	emg_status	49
5.5.4	Variable Documentation	50
5.5.4.1	c_mem	50
5.5.4.2	dev_pwm_limit	50
5.5.4.3	dev_tension	50
5.5.4.4	dev_tension_f	50
5.5.4.5	forced_open	50
5.5.4.6	g_measOld	51
5.5.4.7	g_refOld	51
5.5.4.8	g_rx	51
5.5.4.9	interrupt_flag	51
5.5.4.10	master_mode	51
5.5.4.11	pow_tension	51
5.5.4.12	pwm_sign	51
5.5.4.13	receive_meas_from_hand	51
5.5.4.14	reset_last_value_flag	52
5.5.4.15	rest_enabled	52
5.5.4.16	rest_pos_curr_ref	52
5.5.4.17	tau_feedback	52
5.5.4.18	tension_valid	52
5.5.4.19	timer_value	52
5.5.4.20	timer_value0	52
5.5.4.21	watchdog_flag	52
5.6	interruptions.c File Reference	53
5.6.1	Detailed Description	54
5.6.2	Function Documentation	54
5.6.2.1	analog_read_end()	54
5.6.2.2	command_slave()	54

5.6.2.3	deactivate_slave()	54
5.6.2.4	encoder_reading()	54
5.6.2.5	function_scheduler()	55
5.6.2.6	interrupt_manager()	55
5.6.2.7	motor_control()	55
5.6.2.8	overcurrent_control()	55
5.6.2.9	pwm_limit_search()	55
5.6.3	Variable Documentation	55
5.6.3.1	pwm_preload_values	56
5.7	interruptions.h File Reference	56
5.7.1	Detailed Description	58
5.7.2	Function Documentation	58
5.7.2.1	analog_read_end()	58
5.7.2.2	command_slave()	58
5.7.2.3	CY_ISR_PROTO() [1/2]	58
5.7.2.4	CY_ISR_PROTO() [2/2]	58
5.7.2.5	deactivate_slave()	59
5.7.2.6	encoder_reading()	59
5.7.2.7	function_scheduler()	59
5.7.2.8	interrupt_manager()	59
5.7.2.9	motor_control()	59
5.7.2.10	overcurrent_control()	59
5.7.2.11	pwm_limit_search()	60
5.8	main.c File Reference	60
5.8.1	Detailed Description	61
5.9	utils.c File Reference	61
5.9.1	Detailed Description	62
5.9.2	Function Documentation	63
5.9.2.1	calc_turns_fcn()	63
5.9.2.2	calibration()	63

5.9.2.3	check_enc_data()	63
5.9.2.4	check_rest_position()	64
5.9.2.5	curr_estim()	64
5.9.2.6	filter_acc_1()	64
5.9.2.7	filter_acc_2()	64
5.9.2.8	filter_acc_3()	65
5.9.2.9	filter_ch1()	65
5.9.2.10	filter_ch2()	66
5.9.2.11	filter_curr_diff()	66
5.9.2.12	filter_i1()	66
5.9.2.13	filter_v()	67
5.9.2.14	filter_vel_1()	67
5.9.2.15	filter_vel_2()	67
5.9.2.16	filter_vel_3()	68
5.9.2.17	filter_voltage()	68
5.9.2.18	my_mod()	68
5.9.2.19	my_round()	69
5.10	utils.h File Reference	69
5.10.1	Detailed Description	71
5.10.2	Function Documentation	71
5.10.2.1	calc_turns_fcn()	71
5.10.2.2	calibration()	71
5.10.2.3	check_enc_data()	72
5.10.2.4	check_rest_position()	72
5.10.2.5	curr_estim()	72
5.10.2.6	filter_acc_1()	72
5.10.2.7	filter_acc_2()	73
5.10.2.8	filter_acc_3()	73
5.10.2.9	filter_ch1()	74
5.10.2.10	filter_ch2()	74
5.10.2.11	filter_curr_diff()	74
5.10.2.12	filter_i1()	75
5.10.2.13	filter_v()	75
5.10.2.14	filter_vel_1()	75
5.10.2.15	filter_vel_2()	76
5.10.2.16	filter_vel_3()	76
5.10.2.17	filter_voltage()	77
5.10.2.18	my_mod()	77
5.10.2.19	my_round()	77

Chapter 1

Firmware

This is the firmware of the SoftHand Pro board.

Version

6.1.1

This is the firmware of the SoftHand Pro board in Master configuration. It reads EMG connected to socket and control a motor of an attached SoftHand. Also can read and convert analog measurements connected to the PSoC microcontroller.

Chapter 2

Data Structure Index

2.1 Data Structures

Here are the data structures with brief descriptions:

st_calib	Hand calibration structure	7
st_data	Data sent/received structure	8
st_dev	Device related structure	9
st_meas	10
st_mem	12
st_ref	Motor Reference structure	18

Chapter 3

File Index

3.1 File List

Here is a list of all documented files with brief descriptions:

command_processing.c	
Command processing functions	19
command_processing.h	
Received commands processing functions	26
commands.h	
Definitions for SoftHand commands, parameters and packages	34
device.h	??
globals.c	
Global variables	39
globals.h	
Global definitions and macros are set in this file	43
interruptions.c	
Interruption handling and firmware core functions	53
interruptions.h	
Interruptions header file	56
main.c	
Firmware main file	60
utils.c	
Definition of utility functions	61
utils.h	
Utility functions declaration	69

Chapter 4

Data Structure Documentation

4.1 st_calib Struct Reference

Hand calibration structure.

```
#include <globals.h>
```

Data Fields

- uint8 **enabled**
- uint8 **direction**
- int16 **speed**
- int16 **repetitions**

4.1.1 Detailed Description

Hand calibration structure.

4.1.2 Field Documentation

4.1.2.1 direction

```
uint8 direction
```

Direction of motor winding.

4.1.2.2 enabled

```
uint8 enabled
```

Calibration enabling flag.

4.1.2.3 repetitions

`int16 repetitions`

Number of cycles of hand closing/opening.

4.1.2.4 speed

`int16 speed`

Speed of hand opening/closing.

The documentation for this struct was generated from the following file:

- **globals.h**

4.2 st_data Struct Reference

Data sent/received structure.

```
#include <globals.h>
```

Data Fields

- `uint8` **buffer** [128]
- `int16` **length**
- `int16` **ind**
- `uint8` **ready**

4.2.1 Detailed Description

Data sent/received structure.

4.2.2 Field Documentation

4.2.2.1 buffer

`uint8 buffer[128]`

Data buffer [CMD | DATA | CHECKSUM].

4.2.2.2 ind

```
int16 ind
```

Data buffer index.

4.2.2.3 length

```
int16 length
```

Data buffer length.

4.2.2.4 ready

```
uint8 ready
```

Data buffer flag to see if the data is ready.

The documentation for this struct was generated from the following file:

- **globals.h**

4.3 st_dev Struct Reference

Device related structure.

```
#include <globals.h>
```

Data Fields

- int32 **tension**
- float **tension_conv_factor**
- int8 **pwm_limit**
- uint8 **tension_valid**

4.3.1 Detailed Description

Device related structure.

4.3.2 Field Documentation

4.3.2.1 pwm_limit

```
int8 pwm_limit
```

Limit on pwm value driven to the motor.

4.3.2.2 tension

```
int32 tension
```

Power supply tension.

4.3.2.3 tension_conv_factor

```
float tension_conv_factor
```

Used to calculate input tension.

4.3.2.4 tension_valid

```
uint8 tension_valid
```

Flag checked if the power supply has a valid value.

The documentation for this struct was generated from the following file:

- **globals.h**

4.4 st_meas Struct Reference

Data Fields

- int32 **pos** [**NUM_OF_SENSORS**]
- int32 **curr** [**NUM_OF_MOTORS**]
- int8 **rot** [**NUM_OF_SENSORS**]
- int32 **emg** [**NUM_OF_EMGS**]
- int32 **vel** [**NUM_OF_SENSORS**]
- int32 **acc** [**NUM_OF_SENSORS**]
- int32 **hand_meas**

4.4.1 Field Documentation

4.4.1.1 acc

```
int32 acc[ NUM_OF_SENSORS]
```

Encoder rotational acceleration.

4.4.1.2 curr

```
int32 curr[ NUM_OF_MOTORS]
```

Motor current and current estimation.

4.4.1.3 emg

```
int32 emg[ NUM_OF_EMGS]
```

EMG sensors values.

4.4.1.4 hand_meas

```
int32 hand_meas
```

Position measurement from SH.

4.4.1.5 pos

```
int32 pos[ NUM_OF_SENSORS]
```

Encoder sensor position.

4.4.1.6 rot

```
int8 rot[ NUM_OF_SENSORS]
```

Encoder sensor rotations.

4.4.1.7 vel

```
int32 vel[ NUM_OF_SENSORS]
```

Encoder rotational velocity.

The documentation for this struct was generated from the following file:

- **globals.h**

4.5 st_mem Struct Reference

Data Fields

- uint8 **flag**
- uint8 **id**
- int32 **k_p**
- int32 **k_i**
- int32 **k_d**
- int32 **k_p_c**
- int32 **k_i_c**
- int32 **k_d_c**
- int32 **k_p_dl**
- int32 **k_i_dl**
- int32 **k_d_dl**
- int32 **k_p_c_dl**
- int32 **k_i_c_dl**
- int32 **k_d_c_dl**
- uint8 **activ**
- uint8 **input_mode**
- uint8 **control_mode**
- uint8 **res** [NUM_OF_SENSORS]
- int32 **m_off** [NUM_OF_SENSORS]
- float **m_mult** [NUM_OF_SENSORS]
- uint8 **pos_lim_flag**
- int32 **pos_lim_inf** [NUM_OF_MOTORS]
- int32 **pos_lim_sup** [NUM_OF_MOTORS]
- int32 **max_step_pos**
- int32 **max_step_neg**
- int16 **current_limit**
- uint16 **emg_threshold** [NUM_OF_EMGS]
- uint8 **emg_calibration_flag**
- uint32 **emg_max_value** [NUM_OF_EMGS]
- uint8 **emg_speed**
- uint8 **double_encoder_on_off**
- int8 **motor_handle_ratio**
- uint8 **activate_pwm_rescaling**
- float **curr_lookup** [LOOKUP_DIM]
- uint8 **baud_rate**
- uint8 **watchdog_period**
- int32 **rest_pos**
- float **rest_delay**
- float **rest_vel**
- float **rest_ratio**

4.5.1 Field Documentation

4.5.1.1 `activ`

`uint8 activ`

Startup activation.

4.5.1.2 `activate_pwm_rescaling`

`uint8 activate_pwm_rescaling`

Activation of PWM rescaling for 12V motors.

4.5.1.3 `baud_rate`

`uint8 baud_rate`

Baud Rate set.

4.5.1.4 `control_mode`

`uint8 control_mode`

Motor Control mode.

4.5.1.5 `curr_lookup`

`float curr_lookup[LOOKUP_DIM]`

Table of values to get estimated curr.

4.5.1.6 `current_limit`

`int16 current_limit`

Limit for absorbed current.

4.5.1.7 `double_encoder_on_off`

`uint8 double_encoder_on_off`

Double encoder ON/OFF.

4.5.1.8 `emg_calibration_flag`

`uint8 emg_calibration_flag`

Enable emg calibration on startup.

4.5.1.9 emg_max_value

```
uint32 emg_max_value[ NUM_OF_EMGS]
```

Maximum value for EMG.

4.5.1.10 emg_speed

```
uint8 emg_speed
```

Maximum closure speed when using EMG.

4.5.1.11 emg_threshold

```
uint16 emg_threshold[ NUM_OF_EMGS]
```

Minimum value for activation of EMG control.

4.5.1.12 flag

```
uint8 flag
```

If checked the device has been configured.

4.5.1.13 id

```
uint8 id
```

Device id.

4.5.1.14 input_mode

```
uint8 input_mode
```

Motor Input mode.

4.5.1.15 k_d

```
int32 k_d
```

Position controller derivative constant.

4.5.1.16 k_d_c

```
int32 k_d_c
```

Current controller derivative constant.

4.5.1.17 k_d_c_dl

```
int32 k_d_c_dl
```

Double loop current controller deriv. constant.

4.5.1.18 k_d_dl

```
int32 k_d_dl
```

Double loop position controller deriv. constant.

4.5.1.19 k_i

```
int32 k_i
```

Position controller integrative constant.

4.5.1.20 k_i_c

```
int32 k_i_c
```

Current controller integrative constant.

4.5.1.21 k_i_c_dl

```
int32 k_i_c_dl
```

Double loop current controller integr. constant.

4.5.1.22 k_i_dl

```
int32 k_i_dl
```

Double loop position controller integr. constant.

4.5.1.23 k_p

```
int32 k_p
```

Position controller proportional constant.

4.5.1.24 k_p_c

```
int32 k_p_c
```

Current controller proportional constant.

4.5.1.25 k_p_c_dl

```
int32 k_p_c_dl
```

Double loop current controller prop. constant.

4.5.1.26 k_p_dl

```
int32 k_p_dl
```

Double loop position controller prop. constant.

4.5.1.27 m_mult

```
float m_mult[ NUM_OF_SENSORS]
```

Measurement multiplier.

4.5.1.28 m_off

```
int32 m_off[ NUM_OF_SENSORS]
```

Measurement offset.

4.5.1.29 max_step_neg

```
int32 max_step_neg
```

Maximum number of steps per cycle for negative positions.

4.5.1.30 max_step_pos

```
int32 max_step_pos
```

Maximum number of steps per cycle for positive positions.

4.5.1.31 motor_handle_ratio

```
int8 motor_handle_ratio
```

Discrete multiplier for handle device.

4.5.1.32 pos_lim_flag

```
uint8 pos_lim_flag
```

Position limit active/inactive.

4.5.1.33 pos_lim_inf

```
int32 pos_lim_inf[ NUM_OF_MOTORS]
```

Inferior position limit for motors.

4.5.1.34 pos_lim_sup

```
int32 pos_lim_sup[ NUM_OF_MOTORS]
```

Superior position limit for motors.

4.5.1.35 res

```
uint8 res[ NUM_OF_SENSORS]
```

Angle resolution.

4.5.1.36 rest_delay

```
float rest_delay
```

Hand rest position delay while in EMG mode.

4.5.1.37 rest_pos

```
int32 rest_pos
```

Hand rest position while in EMG mode.

4.5.1.38 rest_ratio

```
float rest_ratio
```

Hand rest ratio between velocity closure and rest position error.

4.5.1.39 rest_vel

```
float rest_vel
```

Hand velocity closure for rest position reaching.

4.5.1.40 watchdog_period

```
uint8 watchdog_period
```

Watchdog period setted, 255 = disable.

The documentation for this struct was generated from the following file:

- **globals.h**

4.6 st_ref Struct Reference

Motor Reference structure.

```
#include <globals.h>
```

Data Fields

- int32 **pos** [**NUM_OF_MOTORS**]
- int32 **curr** [**NUM_OF_MOTORS**]
- int32 **pwm** [**NUM_OF_MOTORS**]
- uint8 **onoff**

4.6.1 Detailed Description

Motor Reference structure.

4.6.2 Field Documentation

4.6.2.1 curr

```
int32 curr[ NUM_OF_MOTORS]
```

Motor current reference.

4.6.2.2 onoff

```
uint8 onoff
```

Motor drivers enable.

4.6.2.3 pos

```
int32 pos[ NUM_OF_MOTORS]
```

Motor position reference.

4.6.2.4 pwm

```
int32 pwm[ NUM_OF_MOTORS]
```

Motor direct pwm control.

The documentation for this struct was generated from the following file:

- **globals.h**

Chapter 5


File Documentation

5.1 command_processing.c File Reference

Command processing functions.

```
#include <command_processing.h>
#include <interruptions.h>
#include <stdio.h>
#include <utils.h>
#include "commands.h"
```

Include dependency graph for command_processing.c:

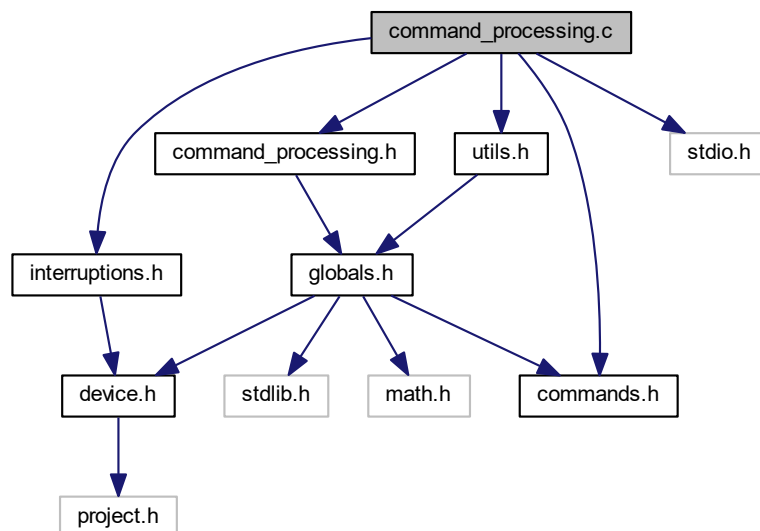

Functions

- void **commProcess** (void)
- void **infoSend** (void)
- void **infoGet** (uint16 info_type)
- void **get_param_list** (uint16 index)
- void **setZeros** ()
- void **infoPrepare** (unsigned char *info_string)
- void **commWrite_old_id** (uint8 *packet_data, uint16 packet_lenght, uint8 old_id)
- void **commWrite** (uint8 *packet_data, uint16 packet_lenght)
- int32 **commReadMeasFromAnother** ()
- void **commWriteAnother** (uint8 *packet_data, uint16 packet_lenght)
- uint8 **LCRChecksum** (uint8 *data_array, uint8 data_length)
- void **sendAcknowledgment** (uint8 value)
- uint8 **memStore** (int displacement)
- void **memRecall** (void)
- uint8 **memRestore** (void)
- uint8 **memInit** (void)
- void **cmd_get_measurements** ()
- void **cmd_get_measurements_from_SH** ()
- void **cmd_set_inputs** ()
- void **cmd_activate** ()
- void **cmd_get_activate** ()
- void **cmd_get_curr_and_meas** ()
- void **cmd_get_currents** ()
- void **cmd_set_baudrate** ()
- void **cmd_ping** ()
- void **cmd_set_watchdog** ()
- void **cmd_get_inputs** ()
- void **cmd_store_params** ()
- void **cmd_get_emg** ()

Variables

- reg8 * **EEPROM_ADDR** = (reg8 *) CYDEV_EE_BASE

5.1.1 Detailed Description

Command processing functions.

Date

October 01, 2017

Author

Centro "E.Piaggio"

Copyright

(C) 2012-2016 qbrobotics. All rights reserved.
 (C) 2017 Centro "E.Piaggio". All rights reserved.

5.1.2 Function Documentation

5.1.2.1 cmd_activate()

```
void cmd_activate ( )
```

This function activates the board

5.1.2.2 cmd_get_activate()

```
void cmd_get_activate ( )
```

This function gets the board activation status and puts it in the package to be sent.

5.1.2.3 cmd_get_currents()

```
void cmd_get_currents ( )
```

This function gets the motor current and puts it in the package to be sent.

5.1.2.4 cmd_get_emg()

```
void cmd_get_emg ( )
```

This function gets the electromyographic sensors measurements and puts them in the package to be sent.

5.1.2.5 cmd_get_inputs()

```
void cmd_get_inputs ( )
```

This function gets the current motor reference inputs and puts them in the package to be sent.

5.1.2.6 cmd_get_measurements()

```
void cmd_get_measurements ( )
```

Bunch of functions used on request from UART communication

5.1.2.7 cmd_get_measurements_from_SH()

```
void cmd_get_measurements_from_SH ( )
```

This function gets the encoders measurements from SH in master mode and puts them in the package to be sent.

5.1.2.8 cmd_ping()

```
void cmd_ping ( )
```

This function is used to ping the device and see if is connected.

5.1.2.9 cmd_set_baudrate()

```
void cmd_set_baudrate ( )
```

This function sets the desired communication baudrate. It is possible to select a value equal to 460800 or 2000000.

5.1.2.10 cmd_set_inputs()

```
void cmd_set_inputs ( )
```

This function gets the inputs from the received package and sets them as motor reference.

5.1.2.11 cmd_set_watchdog()

```
void cmd_set_watchdog ( )
```

This function sets the watchdog timer to the one received from the package. The board automatically deactivate when the time equivalent, to watchdog timer, has passed.

5.1.2.12 cmd_store_params()

```
void cmd_store_params ( )
```

This function stores the parameters to the EEPROM memory

5.1.2.13 commProcess()

```
void commProcess ( )
```

This function unpacks the received package, depending on the command received.

5.1.2.14 commReadMeasFromAnother()

```
int32 commReadMeasFromAnother ( )
```

This function reads on the serial port the measurements package from another board.

5.1.2.15 commWrite()

```
void commWrite (
    uint8 * packet_data,
    uint16 packet_lenght )
```

This function writes on the serial port the package that needs to be sent to the user.

Parameters

<i>packet_data</i>	The array of data that must be written.
<i>packet_lenght</i>	The lenght of the data array.

5.1.2.16 commWrite_old_id()

```
void commWrite_old_id (
    uint8 * packet_data,
    uint16 packet_lenght,
    uint8 old_id )
```

This function writes on the serial port the package that needs to be sent to the user. Is used only when a new is set, to communicate back to the APIs that the new ID setting went fine or there was an error.

Parameters

<i>packet_data</i>	The array of data that must be written.
<i>packet_lenght</i>	The lenght of the data array.
<i>old_id</i>	The previous id of the board, before setting a new one

5.1.2.17 commWriteAnother()

```
void commWriteAnother (
    uint8 * packet_data,
    uint16 packet_lenght )
```

This function writes on the serial port the package that needs to be sent to another board.

Parameters

<i>packet_data</i>	The array of data that must be written.
<i>packet_lenght</i>	The lenght of the data array.

5.1.2.18 get_param_list()

```
void get_param_list (
    uint16 index )
```

This function, depending on the **Firmware** (p. 1) received, gets the list of parameters with their values and sends them to user or sets a parameter from all the parameters of the device.

Parameters

<i>index</i>	The index of the parameters to be setted. If 0 gets full parameters list.
--------------	---

5.1.2.19 infoGet()

```
void infoGet (
    uint16 info_type )
```

This function sends the firmware information prepared with **infoPrepare** (p.32) through the serial port to the user interface. Is used when the ID is specified.

Parameters

<i>info_type</i>	The type of the information needed.
------------------	-------------------------------------

5.1.2.20 infoPrepare()

```
void infoPrepare (
    unsigned char * info_string )
```

This function is used to prepare the information string about the firmware of the device.

Parameters

<i>info_string</i>	An array of chars containing firmware informations.
--------------------	---

5.1.2.21 infoSend()

```
void infoSend ( )
```

This function sends the firmware information prepared with **infoPrepare** (p.32) through the serial port to the user interface. Is used when no ID is specified.

5.1.2.22 LCRChecksum()

```
uint8 LCRChecksum (
    uint8 * data_array,
    uint8 data_length )
```

This function calculates a checksum of the array to see if the received data is consistent.

Parameters

<i>data_array</i>	The array of data that must be checked.
<i>data_lenght</i>	Lenght of the data array that must be checked.

Returns

The calculated checksum for the relative *data_array*.

5.1.2.23 memInit()

```
uint8 memInit ( )
```

This functions initializes the memory. It is used also to restore the the parameters to their default values.

Returns

A true value if the memory is correctly initialized, false otherwise.

5.1.2.24 memRecall()

```
void memRecall ( )
```

This function loads user's settings from the EEPROM.

5.1.2.25 memRestore()

```
uint8 memRestore ( )
```

This function loads default settings from the EEPROM.

Returns

A true value if the memory is correctly restored, false otherwise.

5.1.2.26 memStore()

```
uint8 memStore (
    int displacement )
```

This function stores the setted parameters to the internal EEPROM memory. It is usually called, by the user, after a parameter is set.

Parameters

<i>displacement</i>	The address where the parameters will be written.
---------------------	---

Returns

A true value if the memory is correctly stored, false otherwise.

5.1.2.27 sendAcknowledgment()

```
void sendAcknowledgment (
    uint8 value )
```

This functions sends an acknowledgment to see if a command has been executed properly or not.

Parameters

<i>value</i>	An ACK_OK(1) or ACK_ERROR(0) value.
--------------	-------------------------------------

5.1.2.28 setZeros()

```
void setZeros ( )
```


This function sets the encoders zero position.

5.2 command_processing.h File Reference


Received commands processing functions.

```
#include <globals.h>
```

Include dependency graph for command_processing.h:

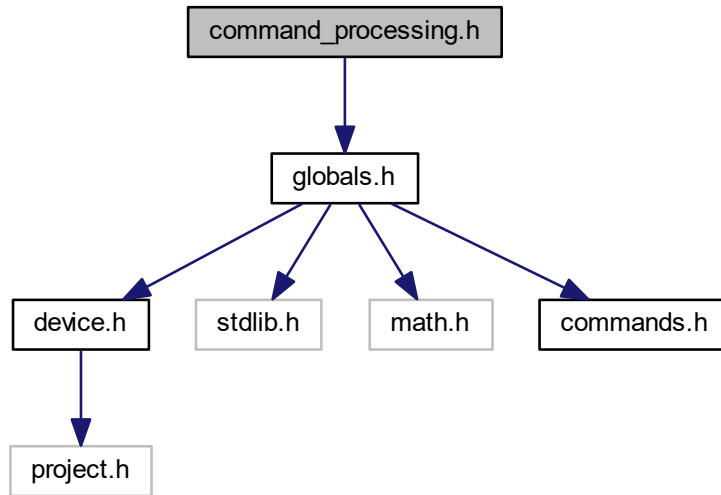

This graph shows which files directly or indirectly include this file:

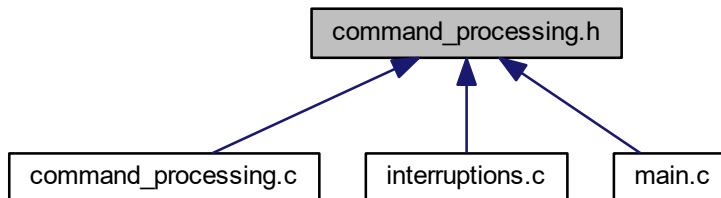

Functions

Firmware information functions

- void **infoPrepare** (unsigned char *info_string)
- void **infoSend** ()
- void **infoGet** (uint16 info_type)

Command receiving and sending functions

- void **commProcess** ()
- void **commWrite_old_id** (uint8 *packet_data, uint16 packet_lenght, uint8 old_id)

- void **commWrite** (uint8 *packet_data, uint16 packet_length)
- void **commWriteAnother** (uint8 *packet_data, uint16 packet_length)
- int32 **commReadMeasFromAnother** ()

Memory management functions

- void **get_param_list** (uint16 index)
- void **setZeros** ()
- uint8 **memStore** (int displacement)
- void **memRecall** ()
- uint8 **memRestore** ()
- uint8 **memInit** ()

Utility functions

- uint8 **LCRCChecksum** (uint8 *data_array, uint8 data_length)
- void **sendAcknowledgment** (uint8 value)

Command processing functions

- void **cmd_activate** ()
- void **cmd_set_inputs** ()
- void **cmd_get_measurements** ()
- void **cmd_get_measurements_from_SH** ()
- void **cmd_get_currents** ()
- void **cmd_get_emg** ()
- void **cmd_set_watchdog** ()
- void **cmd_get_activate** ()
- void **cmd_set_baudrate** ()
- void **cmd_get_inputs** ()
- void **cmd_store_params** ()
- void **cmd_ping** ()

5.2.1 Detailed Description

Received commands processing functions.

Date

October 01, 2017

Author

Centro "E.Piaggio"

Copyright

(C) 2012-2016 qbrobotics. All rights reserved.
(C) 2017 Centro "E.Piaggio". All rights reserved.

This file contains all the definitions of the functions used to process the commands sent from the user interfaces (simulink, command line, GUI)

5.2.2 Function Documentation

5.2.2.1 cmd_activate()

```
void cmd_activate ( )
```

This function activates the board

5.2.2.2 cmd_get_activate()

```
void cmd_get_activate ( )
```

This function gets the board activation status and puts it in the package to be sent.

5.2.2.3 cmd_get_currents()

```
void cmd_get_currents ( )
```

This function gets the motor current and puts it in the package to be sent.

5.2.2.4 cmd_get_emg()

```
void cmd_get_emg ( )
```

This function gets the electromyographic sensors measurements and puts them in the package to be sent.

5.2.2.5 cmd_get_inputs()

```
void cmd_get_inputs ( )
```

This function gets the current motor reference inputs and puts them in the package to be sent.

5.2.2.6 cmd_get_measurements()

```
void cmd_get_measurements ( )
```

This function gets the encoders measurements and puts them in the package to be sent.

Bunch of functions used on request from UART communication

5.2.2.7 cmd_get_measurements_from_SH()

```
void cmd_get_measurements_from_SH ( )
```

This function gets the encoders measurements from SH in master mode and puts them in the package to be sent.

5.2.2.8 cmd_ping()

```
void cmd_ping ( )
```

This function is used to ping the device and see if is connected.

5.2.2.9 cmd_set_baudrate()

```
void cmd_set_baudrate ( )
```

This function sets the desired communication baudrate. It is possible to select a value equal to 460800 or 2000000.

5.2.2.10 cmd_set_inputs()

```
void cmd_set_inputs ( )
```

This function gets the inputs from the received package and sets them as motor reference.

5.2.2.11 cmd_set_watchdog()

```
void cmd_set_watchdog ( )
```

This function sets the watchdog timer to the one received from the package. The board automatically deactivate when the time equivalent, to watchdog timer, has passed.

5.2.2.12 cmd_store_params()

```
void cmd_store_params ( )
```

This function stores the parameters to the EEPROM memory

5.2.2.13 commProcess()

```
void commProcess ( )
```

This function unpacks the received package, depending on the command received.

5.2.2.14 commReadMeasFromAnother()

```
int32 commReadMeasFromAnother ( )
```

This function reads on the serial port the measurements package from another board.

5.2.2.15 commWrite()

```
void commWrite (
    uint8 * packet_data,
    uint16 packet_lenght )
```

This function writes on the serial port the package that needs to be sent to the user.

Parameters

<i>packet_data</i>	The array of data that must be written.
<i>packet_lenght</i>	The lenght of the data array.

5.2.2.16 commWrite_old_id()

```
void commWrite_old_id (
    uint8 * packet_data,
    uint16 packet_lenght,
    uint8 old_id )
```

This function writes on the serial port the package that needs to be sent to the user. Is used only when a new is set, to communicate back to the APIs that the new ID setting went fine or there was an error.

Parameters

<i>packet_data</i>	The array of data that must be written.
<i>packet_lenght</i>	The lenght of the data array.
<i>old_id</i>	The previous id of the board, before setting a new one

5.2.2.17 commWriteAnother()

```
void commWriteAnother (
    uint8 * packet_data,
    uint16 packet_lenght )
```

This function writes on the serial port the package that needs to be sent to another board.

Parameters

<i>packet_data</i>	The array of data that must be written.
<i>packet_lenght</i>	The lenght of the data array.

5.2.2.18 get_param_list()

```
void get_param_list (
    uint16 index )
```

This function, depending on the **Firmware** (p. 1) received, gets the list of parameters with their values and sends them to user or sets a parameter from all the parameters of the device.

Parameters

<i>index</i>	The index of the parameters to be setted. If 0 gets full parameters list.
--------------	---

5.2.2.19 infoGet()

```
void infoGet (
    uint16 info_type )
```

This function sends the firmware information prepared with **infoPrepare** (p.32) through the serial port to the user interface. Is used when the ID is specified.

Parameters

<i>info_type</i>	The type of the information needed.
------------------	-------------------------------------

5.2.2.20 infoPrepare()

```
void infoPrepare (
    unsigned char * info_string )
```

This function is used to prepare the information string about the firmware of the device.

Parameters

<i>info_string</i>	An array of chars containing firmware informations.
--------------------	---

5.2.2.21 infoSend()

```
void infoSend ( )
```

This function sends the firmware information prepared with **infoPrepare** (p.32) through the serial port to the user interface. Is used when no ID is specified.

5.2.2.22 LCRChecksum()

```
uint8 LCRChecksum (
    uint8 * data_array,
    uint8 data_length )
```

This function calculates a checksum of the array to see if the received data is consistent.

Parameters

<i>data_array</i>	The array of data that must be checked.
<i>data_lenght</i>	Lenght of the data array that must be checked.

Returns

The calculated checksum for the relative *data_array*.

5.2.2.23 memInit()

```
uint8 memInit ( )
```

This functions initializes the memory. It is used also to restore the the parameters to their default values.

Returns

A true value if the memory is correctly initialized, false otherwise.

5.2.2.24 memRecall()

```
void memRecall ( )
```

This function loads user's settings from the EEPROM.

5.2.2.25 memRestore()

```
uint8 memRestore ( )
```

This function loads default settings from the EEPROM.

Returns

A true value if the memory is correctly restored, false otherwise.

5.2.2.26 memStore()

```
uint8 memStore (
    int displacement )
```

This function stores the setted parameters to the internal EEPROM memory. It is usually called, by the user, after a parameter is set.

Parameters

<i>displacement</i>	The address where the parameters will be written.
---------------------	---

Returns

A true value if the memory is correctly stored, false otherwise.

5.2.2.27 sendAcknowledgment()

```
void sendAcknowledgment (
    uint8 value )
```

This functions sends an acknowledgment to see if a command has been executed properly or not.

Parameters

<i>value</i>	An ACK_OK(1) or ACK_ERROR(0) value.
--------------	-------------------------------------

5.2.2.28 setZeros()


```
void setZeros ( )
```

This function sets the encoders zero position.

5.3 commands.h File Reference

Definitions for SoftHand commands, parameters and packages.

This graph shows which files directly or indirectly include this file:

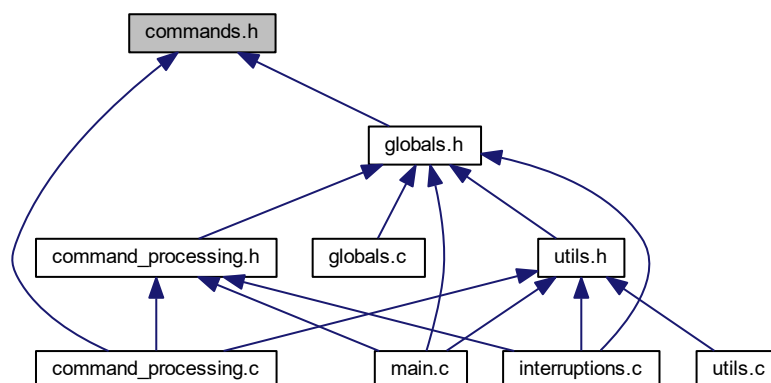

Macros

QB Move Information Strings

- `#define INFO_ALL 0`
All system information.

Enumerations

qbMove and qbHand Commands

- enum `qbmove_command` {
`CMD_PING = 0, CMD_SET_ZEROS = 1, CMD_STORE_PARAMS = 3, CMD_STORE_DEFAULT_P←`
`ARAMS = 4,`
`CMD_RESTORE_PARAMS = 5, CMD_GET_INFO = 6, CMD_SET_VALUE = 7, CMD_GET_VALUE =`
`8,`
`CMD_BOOTLOADER = 9, CMD_INIT_MEM = 10, CMD_CALIBRATE = 11, CMD_GET_PARAM_LIST`
`= 12,`
`CMD_HAND_CALIBRATE = 13, CMD_ACTIVATE = 128, CMD_GET_ACTIVATE = 129, CMD_SET←`
`_INPUTS = 130,`
`CMD_GET_INPUTS = 131, CMD_GET_MEASUREMENTS = 132, CMD_GET_CURRENTS = 133, C←`
`MD_GET_CURR_AND_MEAS = 134,`
`CMD_SET_POS_STIFF = 135, CMD_GET_EMG = 136, CMD_GET_VELOCITIES = 137, CMD_GET←`
`_COUNTERS = 138,`
`CMD_GET_ACCEL = 139, CMD_GET_CURR_DIFF = 140, CMD_SET_CURR_DIFF = 141, CMD_S←`
`ET_CUFF_INPUTS = 142,`
`CMD_SET_WATCHDOG = 143, CMD_SET_BAUDRATE = 144, CMD_EXT_DRIVE = 145, CMD_G←`
`ET_JOYSTICK = 146,`
`CMD_GET_HAND_MEASUREMENTS = 147 }`

qbMove and qbHand Parameters

- `#define PARAM_BYTE_SLOT 50`
Number of bytes reserved to a param information.
- `#define PARAM_MENU_SLOT 150`
Number of bytes reserved to a param menu.
- enum `qbmove_parameter` {
`PARAM_ID = 0, PARAM_PID_CONTROL = 1, PARAM_STARTUP_ACTIVATION = 2, PARAM_INPU←`
`T_MODE = 3,`
`PARAM_CONTROL_MODE = 4, PARAM_MEASUREMENT_OFFSET = 5, PARAM_MEASUREMENT←`
`_MULTIPLIER = 6, PARAM_POS_LIMIT_FLAG = 7,`
`PARAM_POS_LIMIT = 8, PARAM_MAX_STEP_POS = 9, PARAM_MAX_STEP_NEG = 10, PARAM_←`
`POS_RESOLUTION = 11,`
`PARAM_CURRENT_LIMIT = 12, PARAM_EMG_CALIB_FLAG = 13, PARAM_EMG_THRESHOLD = 14,`
`PARAM_EMG_MAX_VALUE = 15,`
`PARAM_EMG_SPEED = 16, PARAM_PID_CURR_CONTROL = 18, PARAM_DOUBLE_ENC_ON_OFF`
`= 19, PARAM_MOT_HANDLE_RATIO = 20,`
`PARAM_MOTOR_SUPPLY = 21, PARAM_CURRENT_LOOKUP = 23, PARAM_DL_POS_PID = 24, P←`
`ARAM_DL_CURR_PID = 25 }`
- enum `qbmove_resolution` {
`RESOLUTION_360 = 0, RESOLUTION_720 = 1, RESOLUTION_1440 = 2, RESOLUTION_2880 = 3,`
`RESOLUTION_5760 = 4, RESOLUTION_11520 = 5, RESOLUTION_23040 = 6, RESOLUTION_46080 = 7,`
`RESOLUTION_92160 = 8 }`

- enum **qbmove_input_mode** {
INPUT_MODE_EXTERNAL = 0, **INPUT_MODE_ENCODER3** = 1, **INPUT_MODE_EMG_PROPORTION**↵
AL = 2, **INPUT_MODE_EMG_INTEGRAL** = 3,
INPUT_MODE_EMG_FCFS = 4, **INPUT_MODE_EMG_FCFS_ADV** = 5 }
- enum **qbmove_control_mode** {
CONTROL_ANGLE = 0, **CONTROL_PWM** = 1, **CONTROL_CURRENT** = 2, **CURR_AND_POS_CONT**↵
ROL = 3,
CONTROL_ANGLE_AND_REST_POS = 4 }
- enum **motor_supply_tipe** { **MAXON_24V** = 0, **MAXON_12V** = 1 }
- enum **acknowledgment_values** { **ACK_ERROR** = 0, **ACK_OK** = 1 }
- enum **data_types** {
TYPE_FLAG = 0, **TYPE_INT8** = 1, **TYPE_UINT8** = 2, **TYPE_INT16** = 3,
TYPE_UINT16 = 4, **TYPE_INT32** = 5, **TYPE_UINT32** = 6, **TYPE_FLOAT** = 7,
TYPE_DOUBLE = 8 }

5.3.1 Detailed Description

Definitions for SoftHand commands, parameters and packages.

Date

October 01, 2017

Author

Centro "E.Piaggio"

Copyright

(C) 2012-2016 qbrobotics. All rights reserved.
 (C) 2017 Centro "E.Piaggio". All rights reserved.

This file is included in the firmware, in its libraries and applications. It contains all definitions that are necessary for the contruction of communication packages.

It includes definitions for all of the device commands, parameters and also the size of answer packages.

5.3.2 Enumeration Type Documentation

Enumerator
5.3.2.1 qbmove_commandenum **qbmove_command****Enumerator**

CMD_PING	Asks for a ping message.
CMD_SET_ZEROS	Command for setting the encoders zero position.
CMD_STORE_PARAMS	Stores all parameters in memory and loads them.
CMD_STORE_DEFAULT_PARAMS	Store current parameters as factory parameters.
CMD_RESTORE_PARAMS	Restore default factory parameters.
CMD_GET_INFO	Asks for a string of information about.
CMD_SET_VALUE	Not Used.
CMD_GET_VALUE	Not Used.
CMD_BOOTLOADER	Sets the bootloader modality to update the firmware.
CMD_INIT_MEM	Initialize the memory with the default values.
CMD_CALIBRATE	Starts the stiffness calibration of the qbMove or the hand closure and opening calibration.
CMD_GET_PARAM_LIST	Command to get the parameters list or to set a defined value chosen by the use.
CMD_HAND_CALIBRATE	Starts a series of opening and closures of the hand.
CMD_ACTIVATE	Command for activating/deactivating the device.
CMD_GET_ACTIVATE	Command for getting device activation state.
CMD_SET_INPUTS	Command for setting reference inputs.
CMD_GET_INPUTS	Command for getting reference inputs.
CMD_GET_MEASUREMENTS	Command for asking device's position measurements.
CMD_GET_CURRENTS	Command for asking device's current measurements.
CMD_GET_CURR_AND_MEAS	Command for asking device's measurements and currents.
CMD_SET_POS_STIFF	Not used in the softHand firmware.
CMD_GET_EMG	Command for asking device's emg sensors measurements.
CMD_GET_VELOCITIES	Command for asking device's velocity measurements.
CMD_GET_COUNTERS	Command for asking device's counters (mostly used for debugging sent commands).
CMD_GET_ACCEL	Command for asking device's acceleration measurements.
CMD_GET_CURR_DIFF	Command for asking device's current difference between a measured one and an estimated one (Only for SoftHand).
CMD_SET_CURR_DIFF	Command used to set current difference modality (Only for Cuff device).
CMD_SET_CUFF_INPUTS	Command used to set Cuff device inputs (Only for Cuff device).
CMD_SET_WATCHDOG	Command for setting watchdog timer or disable it.
CMD_SET_BAUDRATE	Command for setting baudrate of communication.
CMD_EXT_DRIVE	Command to set the actual measurements as inputs to another device (Only for Armslider device).
CMD_GET_JOYSTICK	Command to get the joystick measurements (Only for devices driven by a joystick).
CMD_GET_HAND_MEASUREMENTS	Command for receiving SH position measurements.

5.3.2.2 qbmove_control_mode

enum **qbmove_control_mode**

Enumerator

CONTROL_ANGLE	Classic position control.
CONTROL_PWM	Direct PWM value.
CONTROL_CURRENT	Current control.
CURR_AND_POS_CONTROL	Current and position control.
CONTROL_ANGLE_AND_REST_POS	Position control with rest position check.

5.3.2.3 qbmove_input_mode

enum **qbmove_input_mode**

Enumerator

INPUT_MODE_EXTERNAL	References through external commands (default).
INPUT_MODE_ENCODER3	Encoder 3 drives all inputs.
INPUT_MODE_EMG_PROPORTIONAL	Use EMG measure to proportionally drive the position of the motor 1.
INPUT_MODE_EMG_INTEGRAL	Use 2 EMG signals to drive motor position.
INPUT_MODE_EMG_FCFS	Use 2 EMG. First reaching threshold wins and its value defines hand closure.
INPUT_MODE_EMG_FCFS_ADV	Use 2 EMG. First reaching threshold wins and its value defines hand closure Wait for both EMG to lower under threshold.

5.3.2.4 qbmove_parameter

enum **qbmove_parameter**

Enumerator

PARAM_ID	Device's ID number.
PARAM_PID_CONTROL	PID parameters.
PARAM_STARTUP_ACTIVATION	Start up activation byte.
PARAM_INPUT_MODE	Input mode.
PARAM_CONTROL_MODE	Choose the kind of control between position control, current control, direct PWM value or current+position control.
PARAM_MEASUREMENT_OFFSET	Adds a constant offset to the measurements.
PARAM_MEASUREMENT_MULTIPLIER	Adds a multiplier to the measurements.

Enumerator


PARAM_POS_LIMIT_FLAG	Enable/disable position limiting.
PARAM_POS_LIMIT	Position limit values int32 int32 int32 int32 INF_LIM_1 SUP_LIM_1 INF_LIM_2 SUP_LIM_2
PARAM_MAX_STEP_POS	Used to slow down movements for positive values.
PARAM_MAX_STEP_NEG	Used to slow down movements for negative values.
PARAM_POS_RESOLUTION	Angle resolution for inputs and measurements. Used during communication.
PARAM_CURRENT_LIMIT	Limit for absorbed current.
PARAM_EMG_CALIB_FLAG	Enable calibration on startup.
PARAM_EMG_THRESHOLD	Minimum value to have effect.
PARAM_EMG_MAX_VALUE	Maximum value of EMG.
PARAM_EMG_SPEED	Closure speed when using EMG.
PARAM_PID_CURR_CONTROL	PID current control.
PARAM_DOUBLE_ENC_ON_OFF	Double Encoder Y/N.
PARAM_MOT_HANDLE_RATIO	Multiplier between handle and motor.
PARAM_MOTOR_SUPPLY	Motor supply voltage of the hand.
PARAM_CURRENT_LOOKUP	Table of values used to calculate an estimated current of the SoftHand.
PARAM_DL_POS_PID	Double loop position PID.
PARAM_DL_CURR_PID	Double loop current PID.

5.4 globals.c File Reference

Global variables.

```
#include <globals.h>
```

Include dependency graph for globals.c:

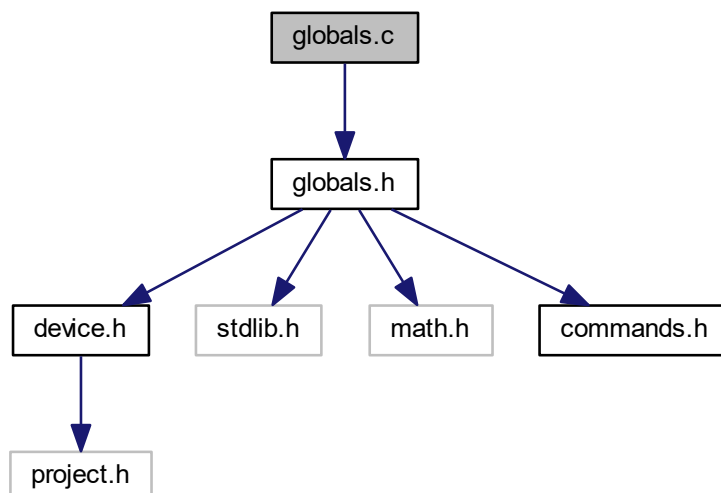

Variables

- struct **st_ref** g_ref g_refNew g_refOld
- struct **st_meas** g_meas g_measOld
- struct **st_data** g_rx
- struct **st_mem** g_mem c_mem
- struct **st_calib** calib
- float **tau_feedback**
- uint32 **timer_value**
- uint32 **timer_value0**
- int32 **dev_tension**
- uint8 **dev_pwm_limit**
- int32 **dev_tension_f**
- int32 **pow_tension**
- CYBIT **reset_last_value_flag**
- CYBIT **tension_valid**
- CYBIT **interrupt_flag** = FALSE
- CYBIT **watchdog_flag**
- int16 **ADC_buf** [4]
- int8 **pwm_sign**
- uint8 **master_mode**
- uint8 **rest_enabled**
- uint8 **forced_open**
- int32 **rest_pos_curr_ref**
- uint8 **receive_meas_from_hand**
- int16 **check2**
- uint8 **check3**
- uint8 **check4**
- int32 **check5**
- int32 **curr_pos_res**

5.4.1 Detailed Description

Global variables.

Date

October 01, 2017

Author

Centro "E.Piaggio"

Copyright

(C) 2012-2016 qbrobotics. All rights reserved.
(C) 2017 Centro "E.Piaggio". All rights reserved.

5.4.2 Variable Documentation

5.4.2.1 c_mem

```
struct st_mem g_mem c_mem
```

Memory parameters.

5.4.2.2 dev_pwm_limit

```
uint8 dev_pwm_limit
```

Device pwm limit. It may change during execution.

5.4.2.3 dev_tension

```
int32 dev_tension
```

Power supply tension.

5.4.2.4 dev_tension_f

```
int32 dev_tension_f
```

Filtered power supply tension.

5.4.2.5 forced_open

```
uint8 forced_open
```

Forced open flag (used in position with rest position control).

5.4.2.6 g_measOld

```
struct st_meas g_meas g_measOld
```

Measurements.

5.4.2.7 g_refOld

```
struct st_ref g_ref g_refNew g_refOld
```

Reference variables.

5.4.2.8 g_rx

```
struct st_data g_rx
```

Incoming/Outcoming data.

5.4.2.9 interrupt_flag

```
CYBIT interrupt_flag = FALSE
```

Interrupt flag enabler.

5.4.2.10 master_mode

```
uint8 master_mode
```

Flag used to set/unset master mode to send messages to other boards.

5.4.2.11 pow_tension

```
int32 pow_tension
```

Computed power supply tension.

5.4.2.12 pwm_sign

```
int8 pwm_sign
```

ADC measurements buffer. Sign of pwm driven. Used to obtain current sign.

5.4.2.13 receive_meas_from_hand

```
uint8 receive_meas_from_hand
```

Position measurement form hand.

5.4.2.14 reset_last_value_flag

```
CYBIT reset_last_value_flag
```

This flag is set when the encoders last values must be resetted.

5.4.2.15 rest_enabled

```
uint8 rest_enabled
```

Rest position flag.

5.4.2.16 rest_pos_curr_ref

```
int32 rest_pos_curr_ref
```

Rest position current reference.

5.4.2.17 tau_feedback

```
float tau_feedback
```

Torque feedback.

5.4.2.18 tension_valid

```
CYBIT tension_valid
```

Tension validation bit.

5.4.2.19 timer_value

```
uint32 timer_value
```

End time of the firmware main loop.

5.4.2.20 timer_value0

```
uint32 timer_value0
```

Start time of the firmware main loop

5.4.2.21 watchdog_flag

```
CYBIT watchdog_flag
```

Watchdog flag enabler.


5.5 globals.h File Reference

Global definitions and macros are set in this file.


```
#include <device.h>
#include "stdlib.h"
#include "math.h"
```

```
#include "commands.h"
```

Include dependency graph for globals.h:

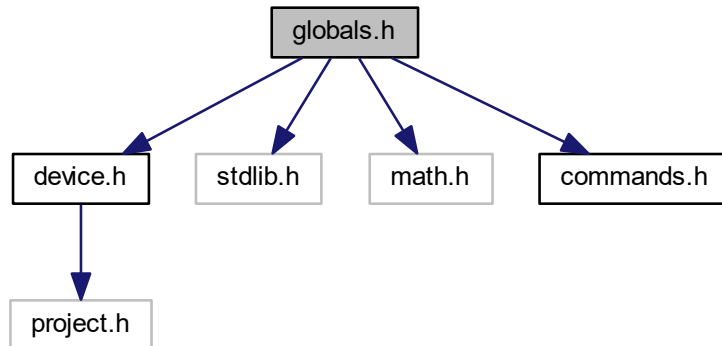

This graph shows which files directly or indirectly include this file:

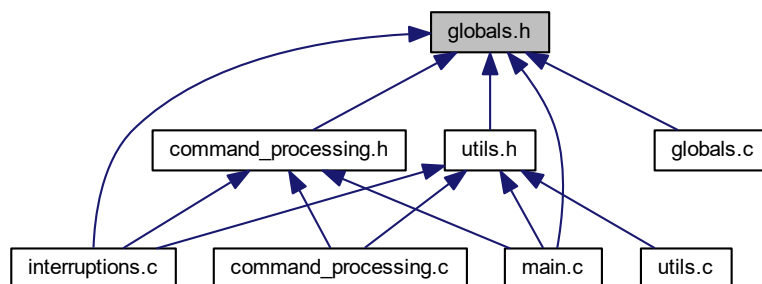

Data Structures

- struct **st_ref**
Motor Reference structure.
- struct **st_meas**
- struct **st_data**
Data sent/received structure.
- struct **st_mem**
- struct **st_dev**
Device related structure.
- struct **st_calib**
Hand calibration structure.

Macros

- `#define VERSION "SH-PRO v6.1.1 - Master w/Rest position"`
- `#define NUM_OF_MOTORS 2`
- `#define NUM_OF_SENSORS 3`
- `#define NUM_OF_EMGS 2`
- `#define NUM_OF_ANALOG_INPUTS 4`
- `#define NUM_OF_PARAMS 25`
- `#define CALIBRATION_DIV 10`
- `#define DIV_INIT_VALUE 1`
- `#define DMA_BYTES_PER_BURST 2`
- `#define DMA_REQUEST_PER_BURST 1`
- `#define DMA_SRC_BASE (CYDEV_PERIPH_BASE)`
- `#define DMA_DST_BASE (CYDEV_SRAM_BASE)`
- `#define WAIT_START 0`
- `#define WAIT_ID 1`
- `#define WAIT_LENGTH 2`
- `#define RECEIVE 3`
- `#define UNLOAD 4`
- `#define FALSE 0`
- `#define TRUE 1`
- `#define DEFAULT_EEPROM_DISPLACEMENT 8`
- `#define MAX_WATCHDOG_TIMER 250`
- `#define PWM_MAX_VALUE 100`
- `#define ANTI_WINDUP 1000`
- `#define DEFAULT_CURRENT_LIMIT 1000`
- `#define CURRENT_HYSTERESIS 10`
- `#define EMG_SAMPLE_TO_DISCARD 500`
- `#define SAMPLES_FOR_MEAN 100`
- `#define SAMPLES_FOR_EMG_MEAN 1000`
- `#define CALIB_DECIMATION 1`
- `#define NUM_OF_CLOSURES 5`
- `#define POS_INTEGRAL_SAT_LIMIT 50000000`
- `#define CURR_INTEGRAL_SAT_LIMIT 100000`
- `#define MIN_CURR_SAT_LIMIT 30`
- `#define LOOKUP_DIM 6`

Enumerations

- `enum emg_status {`
`NORMAL = 0, RESET = 1, DISCARD = 2, SUM_AND_MEAN = 3,`
`WAIT = 4, WAIT_EoC = 5 }`

Variables

- `struct st_ref g_ref g_refNew g_refOld`
- `struct st_meas g_meas g_measOld`
- `struct st_data g_rx`
- `struct st_mem g_mem c_mem`
- `struct st_calib calib`
- `uint32 timer_value`
- `uint32 timer_value0`
- `int32 dev_tension`

- uint8 **dev_pwm_limit**
- int32 **dev_tension_f**
- int32 **pow_tension**
- CYBIT **reset_last_value_flag**
- CYBIT **tension_valid**
- CYBIT **interrupt_flag**
- CYBIT **watchdog_flag**
- float **tau_feedback**
- int16 **ADC_buf** [4]
- int8 **pwm_sign**
- uint8 **master_mode**
- uint8 **rest_enabled**
- uint8 **forced_open**
- int32 **rest_pos_curr_ref**
- uint8 **receive_meas_from_hand**
- int16 **check2**
- uint8 **check3**
- uint8 **check4**
- int32 **check5**
- int32 **curr_pos_res**

5.5.1 Detailed Description

Global definitions and macros are set in this file.

Date

October 01, 2017

Author

Centro "E.Piaggio"

Copyright

(C) 2012-2016 qbrobotics. All rights reserved.
(C) 2017 Centro "E.Piaggio". All rights reserved.

5.5.2 Macro Definition Documentation

5.5.2.1 ANTI_WINDUP

```
#define ANTI_WINDUP 1000
```

Anti windup saturation.

5.5.2.2 CALIBRATION_DIV

```
#define CALIBRATION_DIV 10
```

Frequency divisor for hand calibration (100Hz).

5.5.2.3 CURR_INTEGRAL_SAT_LIMIT

```
#define CURR_INTEGRAL_SAT_LIMIT 100000
```

Anti windup on current control.

5.5.2.4 CURRENT_HYSTERESIS

```
#define CURRENT_HYSTERESIS 10
```

milliAmperes of hysteresis for current control.

5.5.2.5 DEFAULT_CURRENT_LIMIT

```
#define DEFAULT_CURRENT_LIMIT 1000
```

Default Current limit, 0 stands for unlimited.

5.5.2.6 DEFAULT_EEPROM_DISPLACEMENT

```
#define DEFAULT_EEPROM_DISPLACEMENT 8
```

Number of pages occupied by the EEPROM.

5.5.2.7 DIV_INIT_VALUE

```
#define DIV_INIT_VALUE 1
```

Initial value for hand counter calibration.

5.5.2.8 EMG_SAMPLE_TO_DISCARD

```
#define EMG_SAMPLE_TO_DISCARD 500
```

Number of sample to discard before calibration.

5.5.2.9 LOOKUP_DIM

```
#define LOOKUP_DIM 6
```

Dimension of the current lookup table.

5.5.2.10 MAX_WATCHDOG_TIMER

```
#define MAX_WATCHDOG_TIMER 250
```

num * 2 [cs]

5.5.2.11 NUM_OF_ANALOG_INPUTS

```
#define NUM_OF_ANALOG_INPUTS 4
```

Total number of analogic inputs.

5.5.2.12 NUM_OF_EMGS

```
#define NUM_OF_EMGS 2
```

Number of emg channels.

5.5.2.13 NUM_OF_MOTORS

```
#define NUM_OF_MOTORS 2
```

Number of motors.

5.5.2.14 NUM_OF_PARAMS

```
#define NUM_OF_PARAMS 25
```

Number of parameters saved in the EEPROM.

5.5.2.15 NUM_OF_SENSORS

```
#define NUM_OF_SENSORS 3
```

Number of encoders.

5.5.2.16 POS_INTEGRAL_SAT_LIMIT

```
#define POS_INTEGRAL_SAT_LIMIT 50000000
```

Anti windup on position control.

5.5.2.17 PWM_MAX_VALUE

```
#define PWM_MAX_VALUE 100
```

Maximum value of the PWM signal.

5.5.2.18 RECEIVE

```
#define RECEIVE 3
```

Package data receiving status.

5.5.2.19 SAMPLES_FOR_EMG_MEAN

```
#define SAMPLES_FOR_EMG_MEAN 1000
```

Number of samples used to mean emg values.

5.5.2.20 SAMPLES_FOR_MEAN

```
#define SAMPLES_FOR_MEAN 100
```

Number of samples used to mean current values.

5.5.2.21 UNLOAD

```
#define UNLOAD 4
```

Package data flush status.

5.5.2.22 WAIT_ID

```
#define WAIT_ID 1
```

Package ID waiting status.

5.5.2.23 WAIT_LENGTH

```
#define WAIT_LENGTH 2
```

Package lenght waiting status.

5.5.2.24 WAIT_START

```
#define WAIT_START 0
```

Package start waiting status.

5.5.3 Enumeration Type Documentation

5.5.3.1 emg_status

```
enum emg_status
```

Enumerator

NORMAL	Normal execution.
RESET	Reset analog measurements.
DISCARD	Discard first samples to obtain a correct value.
SUM_AND_MEAN	Sum and mean a definite value of samples.
WAIT	The second emg waits until the first emg has a valid value.
WAIT_EoC	The second emg waits for end of calibration.

5.5.4 Variable Documentation**5.5.4.1 c_mem**

```
struct st_mem g_mem c_mem
```

Memory parameters.

5.5.4.2 dev_pwm_limit

```
uint8 dev_pwm_limit
```

Device pwm limit.

Device pwm limit. It may change during execution.

5.5.4.3 dev_tension

```
int32 dev_tension
```

Power supply tension.

5.5.4.4 dev_tension_f

```
int32 dev_tension_f
```

Filtered power supply tension.

5.5.4.5 forced_open

```
uint8 forced_open
```

Forced open flag (used in position with rest position control).

5.5.4.6 g_measOld

```
struct st_meas g_meas g_measOld
```

Measurements.

5.5.4.7 g_refOld

```
struct st_ref g_ref g_refNew g_refOld
```

Reference variables.

5.5.4.8 g_rx

```
struct st_data g_rx
```

Incoming/Outcoming data.

5.5.4.9 interrupt_flag

```
CYBIT interrupt_flag
```

Interrupt flag enabler.

5.5.4.10 master_mode

```
uint8 master_mode
```

Flag used to set/unset master mode to send messages to other boards.

5.5.4.11 pow_tension

```
int32 pow_tension
```

Computed power supply tension.

5.5.4.12 pwm_sign

```
int8 pwm_sign
```

ADC measurements buffer. Sign of pwm driven. Used to obtain current sign.

5.5.4.13 receive_meas_from_hand

```
uint8 receive_meas_from_hand
```

Position measurement form hand.

5.5.4.14 reset_last_value_flag

```
CYBIT reset_last_value_flag
```

This flag is set when the encoders last values must be resetted.

5.5.4.15 rest_enabled

```
uint8 rest_enabled
```

Rest position flag.

5.5.4.16 rest_pos_curr_ref

```
int32 rest_pos_curr_ref
```

Rest position current reference.

5.5.4.17 tau_feedback

```
float tau_feedback
```

Torque feedback.

5.5.4.18 tension_valid

```
CYBIT tension_valid
```

Tension validation bit.

5.5.4.19 timer_value

```
uint32 timer_value
```

End time of the firmware main loop.

5.5.4.20 timer_value0

```
uint32 timer_value0
```

Start time of the firmware main loop

5.5.4.21 watchdog_flag

```
CYBIT watchdog_flag
```


Watchdog flag enabler.

5.6 interruptions.c File Reference

Interrupt handling and firmware core functions.

```
#include <interruptions.h>
#include <command_processing.h>
#include <globals.h>
#include <utils.h>
```

Include dependency graph for interruptions.c:

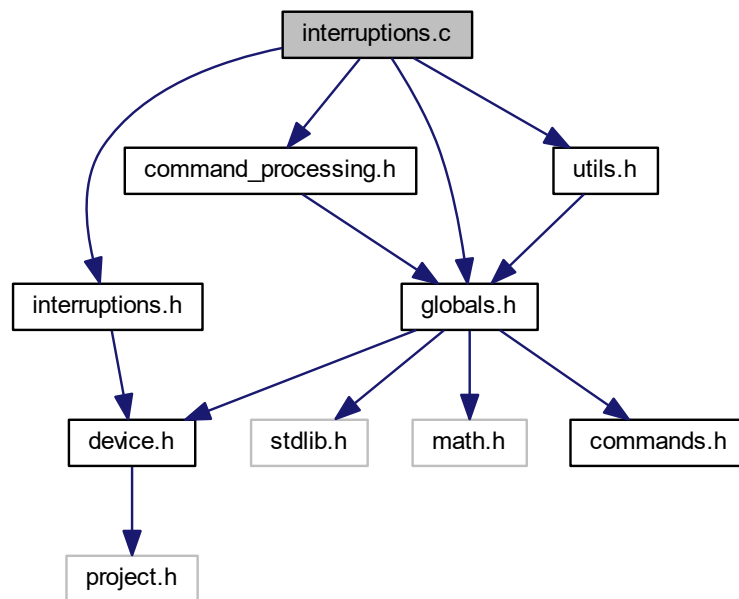

Functions

- **CY_ISR** (ISR_WATCHDOG_Handler)
- **CY_ISR** (ISR_RS485_RX_ExInterrupt)
- void **interrupt_manager** ()
- void **function_scheduler** (void)
- void **motor_control** ()
- void **encoder_reading** (const uint8 idx)
- void **analog_read_end** ()
- void **overcurrent_control** ()
- void **pwm_limit_search** ()
- void **command_slave** ()
- void **deactivate_slave** ()

Variables

- static const uint8 **pwm_preload_values** [29]

5.6.1 Detailed Description

Interruption handling and firmware core functions.

Date

October 01, 2017

Author

Centro "E.Piaggio"

Copyright

(C) 2012-2016 qbrobotics. All rights reserved.

(C) 2017 Centro "E.Piaggio". All rights reserved.

5.6.2 Function Documentation

5.6.2.1 analog_read_end()

```
void analog_read_end ( )
```

This function executes and terminates the analog readings.

5.6.2.2 command_slave()

```
void command_slave ( )
```

This function is used to create a package and send it to another device, only if the actual board is in master mode.

5.6.2.3 deactivate_slave()

```
void deactivate_slave ( )
```

This function is used to create a package to deactivate motors on another device, only if the actual board is in master mode.

5.6.2.4 encoder_reading()

```
void encoder_reading (
    const uint8 index )
```

This functions reads the value from the encoder pointed by index.

Parameters

<i>index</i>	The number of the encoder that must be read.
--------------	--

5.6.2.5 function_scheduler()

```
void function_scheduler (
    void )
```

This function schedules the other functions in an order that optimizes the controller usage.

5.6.2.6 interrupt_manager()

```
void interrupt_manager ( )
```

This function is called in predefined moments during firmware execution in order to unpack the received package.

5.6.2.7 motor_control()

```
void motor_control ( )
```

This function controls the motor direction and velocity, depending on the input and control modality set.

5.6.2.8 overcurrent_control()

```
void overcurrent_control ( )
```

This function increases or decreases the pwm value, depending on the current absorbed by the motor.

5.6.2.9 pwm_limit_search()

```
void pwm_limit_search ( )
```

This function scales the pwm value of the motor, depending on the power supply voltage, in order to not make the motor wind too fast.

5.6.3 Variable Documentation

5.6.3.1 pwm_preload_values

```
const uint8 pwm_preload_values[29] [static]
```

Initial value:


```
= {100,  
83,  
78,  
76,  
74,  
72,  
70,  
68,  
67,  
65,  
64,  
63,  
62,  
61,  
60,  
59,  
58,  
57,  
56,  
56,  
55,  
54,  
54,  
53,  
52,  
52,  
52,  
51,  
51}
```

5.7 interruptions.h File Reference


Interruptions header file.

```
#include <device.h>
```

Include dependency graph for interruptions.h:

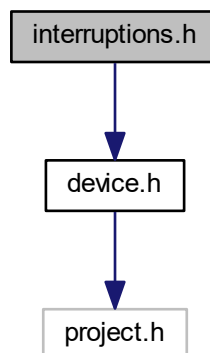

This graph shows which files directly or indirectly include this file:

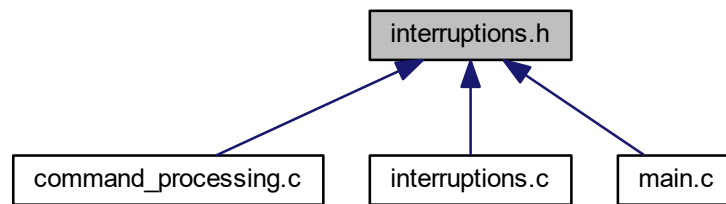

Functions

Interruptions

- **CY_ISR_PROTO** (ISR_RS485_RX_ExInterrupt)
- **CY_ISR_PROTO** (ISR_WATCHDOG_Handler)

General function scheduler

- void **function_scheduler** (void)

Encoder reading function

- void **encoder_reading** (const uint8 index)

Motor control function

- void **motor_control** ()

Analog readings

- void **analog_read_end** ()

Interrupt manager

- void **interrupt_manager** ()

Utility functions

- void **pwm_limit_search** ()
- void **overcurrent_control** ()
- void **command_slave** ()
- void **deactivate_slave** ()

5.7.1 Detailed Description

Interruptions header file.

Date

October 01, 2017

Author

Centro "E.Piaggio"

Copyright

(C) 2012-2016 qbrobotics. All rights reserved.

(C) 2017 Centro "E.Piaggio". All rights reserved.

5.7.2 Function Documentation

5.7.2.1 analog_read_end()

```
void analog_read_end ( )
```

This function executes and terminates the analog readings.

5.7.2.2 command_slave()

```
void command_slave ( )
```

This function is used to create a package and send it to another device, only if the actual board is in master mode.

5.7.2.3 CY_ISR_PROTO() [1/2]

```
CY_ISR_PROTO (
    ISR_RS485_RX_ExInterrupt )
```

This interruption sets a flag to let the firmware know that a communication interruption is pending and needs to be handled. The interruption will be handled in predefined moments during the firmware execution. When this interruption is handled, it unpacks the package received on the RS485 communication bus.

5.7.2.4 CY_ISR_PROTO() [2/2]

```
CY_ISR_PROTO (
    ISR_WATCHDOG_Handler )
```

This interruption sets a flag to let the firmware know that a watchdog interruption is pending and needs to be handled. The interruption will be handled in predefined moments during the firmware execution. When this interruption is handled, it deactivates the board because the watchdog timer has expired.

5.7.2.5 deactivate_slave()

```
void deactivate_slave ( )
```

This function is used to create a package to deactivate motors on another device, only if the actual board is in master mode.

5.7.2.6 encoder_reading()

```
void encoder_reading (
    const uint8 index )
```

This functions reads the value from the encoder pointed by index.

Parameters

<i>index</i>	The number of the encoder that must be read.
--------------	--

5.7.2.7 function_scheduler()

```
void function_scheduler (
    void )
```

This function schedules the other functions in an order that optimizes the controller usage.

5.7.2.8 interrupt_manager()

```
void interrupt_manager ( )
```

This function is called in predefined moments during firmware execution in order to unpack the received package.

5.7.2.9 motor_control()

```
void motor_control ( )
```

This function controls the motor direction and velocity, depending on the input and control modality set.

5.7.2.10 overcurrent_control()

```
void overcurrent_control ( )
```

This function increases or decreases the pwm value, depending on the current absorbed by the motor.

5.7.2.11 pwm_limit_search()

```
void pwm_limit_search ( )
```


This function scales the pwm value of the motor, depending on the power supply voltage, in order to not make the motor wind too fast.

5.8 main.c File Reference

Firmware main file.

```
#include <device.h>
#include <globals.h>
#include <interruptions.h>
#include <command_processing.h>
#include <utils.h>
```

Include dependency graph for main.c:

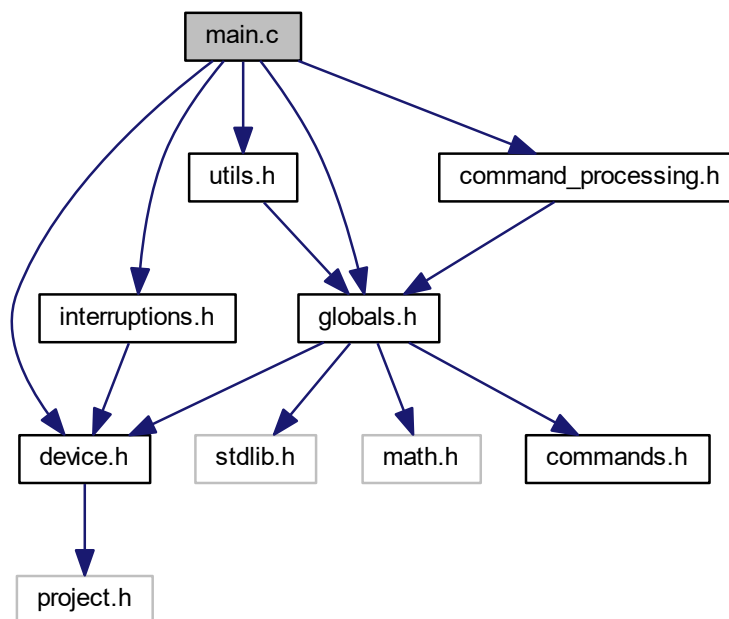

Functions

- `int main ()`

5.8.1 Detailed Description

Firmware main file.

Date

October 01, 2017

Author

Centro "E.Piaggio"

Copyright

(C) 2012-2016 qbrobotics. All rights reserved.

(C) 2017 Centro "E.Piaggio". All rights reserved.


5.9 utils.c File Reference

Definition of utility functions.

```
#include <utils.h>
```

```
#include <math.h>
```

Include dependency graph for utils.c:

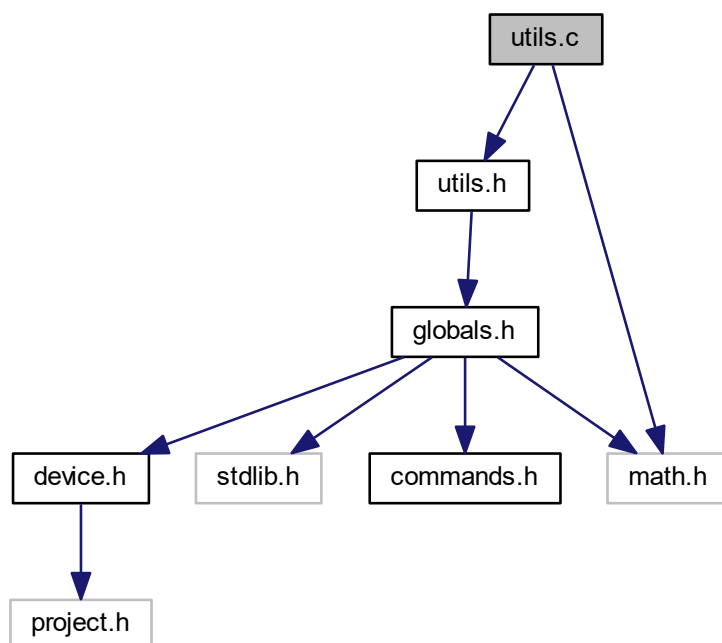

Macros

- **#define N1 15**
Teeth of the first encoder wheel.
- **#define N2 14**
Teeth of the second encoder wheel.
- **#define I1 1**
First wheel invariant value.
- **#define I2 (-1)**
Second wheel invariant value.
- **#define M 65536**
Number of encoder ticks per turn.

Functions

- int32 **curr_estim** (int32 pos, int32 vel, int32 acc)
- int32 **filter_v** (int32 new_value)
- int32 **filter_i1** (int32 new_value)
- int32 **filter_ch1** (int32 new_value)
- int32 **filter_ch2** (int32 new_value)
- int32 **filter_curr_diff** (int32 curr_diff)
- int32 **filter_vel_1** (int32 new_value)
- int32 **filter_vel_2** (int32 new_value)
- int32 **filter_vel_3** (int32 new_value)
- int32 **filter_acc_1** (int32 new_value)
- int32 **filter_acc_2** (int32 new_value)
- int32 **filter_acc_3** (int32 new_value)
- int32 **filter_voltage** (int32 new_value)
- CYBIT **check_enc_data** (const uint32 *value)
- int **my_round** (const double x)
- uint32 **my_mod** (int32 val, int32 divisor)
- void **calibration** (void)
- int **calc_turns_fcn** (const int32 pos1, const int32 pos2)
- void **check_rest_position** (void)

5.9.1 Detailed Description

Definition of utility functions.

Date

October 01, 2017

Author

Centro "E.Piaggio"

Copyright

(C) 2012-2016 qrobotics. All rights reserved.
(C) 2017 Centro "E.Piaggio". All rights reserved.

5.9.2 Function Documentation

5.9.2.1 calc_turns_fcn()

```
int calc_turns_fcn (
    const int32 pos1,
    const int32 pos2 )
```

This function is used at startup to reconstruct the correct turn of the shaft connected to the motor. It need two encoders to work.

Parameters

<i>pos1</i>	First encoder position
<i>pos2</i>	Second encoder position

Returns

Returns the number of turns of motor pulley at startup

5.9.2.2 calibration()

```
void calibration ( )
```

This function counts a series of hand opening and closing used to execute a calibration of the device.

5.9.2.3 check_enc_data()

```
CYBIT check_enc_data (
    const uint32 * value )
```

This function controls if the read encoder data is correct or not.

Parameters

<i>value</i>	A pointer to the encoder data read
--------------	------------------------------------

Returns

Returns 1 if the read data is correct, 0 otherwise

5.9.2.4 check_rest_position()

```
void check_rest_position ( )
```

This function checks for rest position and, in case, gives a position reference to the hand.

5.9.2.5 curr_estim()

```
int32 curr_estim (
    int32 pos,
    int32 vel,
    int32 accel )
```

Function used to obtain current estimation through current lookup table.

Parameters

<i>pos</i>	Position of the encoder in ticks.
<i>vel</i>	Speed of the encoder.
<i>accel</i>	Acceleration of the encoder

Returns

Returns an estimation of the motor current, depending on its position, velocity and acceleration.

5.9.2.6 filter_acc_1()

```
int32 filter_acc_1 (
    int32 value )
```

Filter on first encoder rotational acceleration. The weighted average between the old value and the new one is executed.

Parameters

<i>value</i>	New value of the filter.
--------------	--------------------------

Returns

Returns the filtered first encoder rotational acceleration value

5.9.2.7 filter_acc_2()

```
int32 filter_acc_2 (
    int32 value )
```

Filter on second encoder rotation acceleration (if present). The weighted average between the old value and the new one is executed.

Parameters

<i>value</i>	New value of the filter.
--------------	--------------------------

Returns

Returns the filtered second encoder rotational acceleration value

5.9.2.8 filter_acc_3()

```
int32 filter_acc_3 (  
    int32 value )
```

Filter on third encoder rotation acceleration (if present). The weighted average between the old value and the new one is executed.

Parameters

<i>value</i>	New value of the filter.
--------------	--------------------------

Returns

Returns the filtered third encoder rotational acceleration value

5.9.2.9 filter_ch1()

```
int32 filter_ch1 (  
    int32 value )
```

Filter on the first EMG sensor converted value. The weighted average between the old value and the new one is executed.

Parameters

<i>value</i>	New value of the filter.
--------------	--------------------------

Returns

Returns the filtered emg sensor value

5.9.2.10 filter_ch2()

```
int32 filter_ch2 (  
    int32 value )
```

Filter on the second EMG sensor converted value. The weighted average between the old value and the new one is executed.

Parameters

<i>value</i>	New value of the filter.
--------------	--------------------------

Returns

Returns the filtered emg sensor value

5.9.2.11 filter_curr_diff()

```
int32 filter_curr_diff (  
    int32 curr_diff )
```

Low pass filter on current difference between measured and estimated current

Parameters

<i>curr_diff</i>	Difference between the measured current and the estimated one.
------------------	--

Returns

Returns the filtered current difference value

5.9.2.12 filter_i1()

```
int32 filter_i1 (  
    int32 value )
```

Filter on the motor current converted value. The weighted average between the old value and the new one is executed.

Parameters

<i>value</i>	New value of the filter.
--------------	--------------------------

Returns

Returns the filtered current value

5.9.2.13 filter_v()

```
int32 filter_v (
    int32 new_value )
```

Filter on the converted voltage value. The weighted average between the old value and the new one is executed.

Parameters

<i>new_value</i>	New value of the filter.
------------------	--------------------------

Returns

Returns the filtered voltage value

5.9.2.14 filter_vel_1()

```
int32 filter_vel_1 (
    int32 value )
```

Filter on first encoder rotational speed. The weighted average between the old value and the new one is executed.

Parameters

<i>value</i>	New value of the filter.
--------------	--------------------------

Returns

Returns the filtered first encoder rotational speed value

5.9.2.15 filter_vel_2()

```
int32 filter_vel_2 (
    int32 value )
```

Filter on second encoder rotational speed (if present). The weighted average between the old value and the new one is executed.

Parameters

<i>value</i>	New value of the filter.
--------------	--------------------------

Returns

Returns the filtered second encoder rotational speed value

5.9.2.16 filter_vel_3()

```
int32 filter_vel_3 (  
    int32 value )
```

Filter on third encoder rotational speed (if present). The weighted average between the old value and the new one is executed.

Parameters

<i>value</i>	New value of the filter.
--------------	--------------------------

Returns

Returns the filtered third encoder rotational speed value

5.9.2.17 filter_voltage()

```
int32 filter_voltage (  
    int32 value )
```

Filter on voltage readings. The weighted average between the old value and the new one is executed.

Parameters

<i>value</i>	New value of the filter.
--------------	--------------------------

Returns

Returns the filtered voltage value

5.9.2.18 my_mod()

```
uint32 my_mod (  
    int32 val,  
    int32 divisor )
```


This function computes the module function, returning positive values regardless of whether the value passed is negative

Parameters

<i>val</i>	The value of which the module needs to be calculated
<i>divisor</i>	The divisor according to which the module is calculated

5.9.2.19 my_round()

```
int my_round (
    const double x )
```

This function approximates the value passed to the nearest integer


Parameters

<i>x</i>	The floating point value that needs to be rounded
----------	---


5.10 utils.h File Reference

Utility functions declaration.

```
#include <globals.h>
Include dependency graph for utils.h:
```

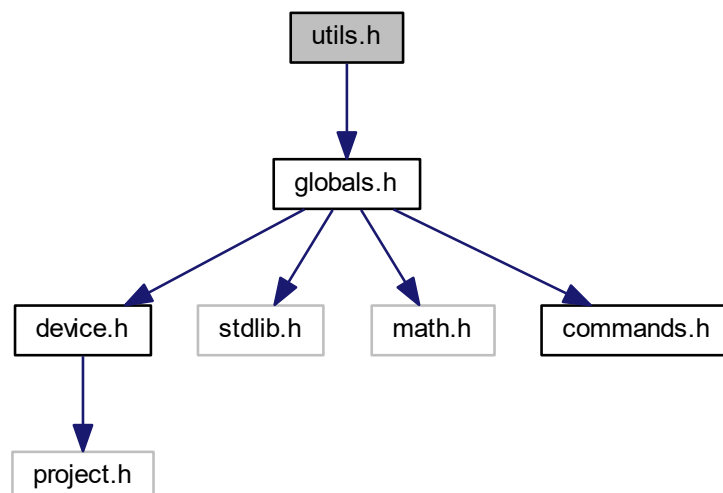

This graph shows which files directly or indirectly include this file:

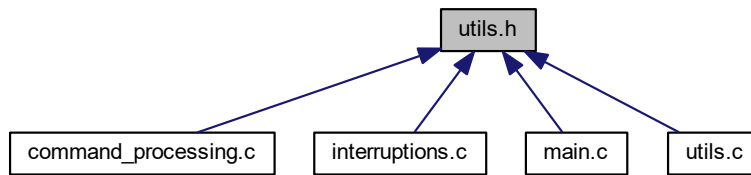

Macros

- `#define TIMER_CLOCK 10000`
- `#define ALPHA 32`
Voltage and current filters constant.
- `#define BETA 50`
Emg filters constant.
- `#define GAMMA 32`
Velocity filters constant.
- `#define DELTA 32`
Acceleration filters constant.
- `#define EPSILON 8`
Voltage readings filter.
- `#define SIGN(A) (((A) >=0) ? (1) : (-1))`
Sign calculation function.

Functions

- `int32 filter_voltage (int32 value)`

Filters

- `int32 filter_v (int32 new_value)`
- `int32 filter_ch1 (int32 value)`
- `int32 filter_ch2 (int32 value)`
- `int32 filter_i1 (int32 value)`
- `int32 filter_vel_1 (int32 value)`
- `int32 filter_vel_2 (int32 value)`
- `int32 filter_vel_3 (int32 value)`
- `int32 filter_acc_1 (int32 value)`
- `int32 filter_acc_2 (int32 value)`
- `int32 filter_acc_3 (int32 value)`

Estimating current and difference

- `int32 curr_estim (int32 pos, int32 vel, int32 accel)`
- `int32 filter_curr_diff (int32 curr_diff)`

Utility functions

- `int my_round (const double x)`
- `uint32 my_mod (int32 val, int32 divisor)`
- `CYBIT check_enc_data (const uint32 *value)`
- `int calc_turns_fcn (const int32 pos1, const int32 pos2)`
- `void calibration ()`
- `void check_rest_position ()`

5.10.1 Detailed Description

Utility functions declaration.

Date

October 01, 2017

Author

Centro "E.Piaggio"

Copyright

(C) 2012-2016 qbrobotics. All rights reserved.

(C) 2017 Centro "E.Piaggio". All rights reserved.

5.10.2 Function Documentation

5.10.2.1 calc_turns_fcn()

```
int calc_turns_fcn (  
    const int32 pos1,  
    const int32 pos2 )
```

This function is used at startup to reconstruct the correct turn of the shaft connected to the motor. It need two encoders to work.

Parameters

<i>pos1</i>	First encoder position
<i>pos2</i>	Second encoder position

Returns

Returns the number of turns of motor pulley at startup

5.10.2.2 calibration()

```
void calibration ( )
```

This function counts a series of hand opening and closing used to execute a calibration of the device.

5.10.2.3 check_enc_data()

```
CYBIT check_enc_data (
    const uint32 * value )
```

This function controls if the read encoder data is correct or not.

Parameters

<i>value</i>	A pointer to the encoder data read
--------------	------------------------------------

Returns

Returns 1 if the read data is correct, 0 otherwise

5.10.2.4 check_rest_position()

```
void check_rest_position ( )
```

This function checks for rest position and, in case, gives a position reference to the hand.

5.10.2.5 curr_estim()

```
int32 curr_estim (
    int32 pos,
    int32 vel,
    int32 accel )
```

Function used to obtain current estimation through current lookup table.

Parameters

<i>pos</i>	Position of the encoder in ticks.
<i>vel</i>	Speed of the encoder.
<i>accel</i>	Acceleration of the encoder

Returns

Returns an estimation of the motor current, depending on its position, velocity and acceleration.

5.10.2.6 filter_acc_1()

```
int32 filter_acc_1 (
    int32 value )
```

Filter on first encoder rotational acceleration. The weighted average between the old value and the new one is executed.

Parameters

<i>value</i>	New value of the filter.
--------------	--------------------------

Returns

Returns the filtered first encoder rotational acceleration value

5.10.2.7 filter_acc_2()

```
int32 filter_acc_2 (  
    int32 value )
```

Filter on second encoder rotation acceleration (if present). The weighted average between the old value and the new one is executed.

Parameters

<i>value</i>	New value of the filter.
--------------	--------------------------

Returns

Returns the filtered second encoder rotational acceleration value

5.10.2.8 filter_acc_3()

```
int32 filter_acc_3 (  
    int32 value )
```

Filter on third encoder rotation acceleration (if present). The weighted average between the old value and the new one is executed.

Parameters

<i>value</i>	New value of the filter.
--------------	--------------------------

Returns

Returns the filtered third encoder rotational acceleration value

5.10.2.9 filter_ch1()

```
int32 filter_ch1 (  
    int32 value )
```

Filter on the first EMG sensor converted value. The weighted average between the old value and the new one is executed.

Parameters

<i>value</i>	New value of the filter.
--------------	--------------------------

Returns

Returns the filtered emg sensor value

5.10.2.10 filter_ch2()

```
int32 filter_ch2 (  
    int32 value )
```

Filter on the second EMG sensor converted value. The weighted average between the old value and the new one is executed.

Parameters

<i>value</i>	New value of the filter.
--------------	--------------------------

Returns

Returns the filtered emg sensor value

5.10.2.11 filter_curr_diff()

```
int32 filter_curr_diff (  
    int32 curr_diff )
```

Low pass filter on current difference between measured and estimated current

Parameters

<i>curr_diff</i>	Difference between the measured current and the estimated one.
------------------	--

Returns

Returns the filtered current difference value

5.10.2.12 filter_i1()

```
int32 filter_i1 (
    int32 value )
```

Filter on the motor current converted value. The weighted average between the old value and the new one is executed.

Parameters

<i>value</i>	New value of the filter.
--------------	--------------------------

Returns

Returns the filtered current value

5.10.2.13 filter_v()

```
int32 filter_v (
    int32 new_value )
```

Filter on the converted voltage value. The weighted average between the old value and the new one is executed.

Parameters

<i>new_value</i>	New value of the filter.
------------------	--------------------------

Returns

Returns the filtered voltage value

5.10.2.14 filter_vel_1()

```
int32 filter_vel_1 (
    int32 value )
```

Filter on first encoder rotational speed. The weighted average between the old value and the new one is executed.

Parameters

<i>value</i>	New value of the filter.
--------------	--------------------------

Returns

Returns the filtered first encoder rotational speed value

5.10.2.15 filter_vel_2()

```
int32 filter_vel_2 (  
    int32 value )
```

Filter on second encoder rotational speed (if present). The weighted average between the old value and the new one is executed.

Parameters

<i>value</i>	New value of the filter.
--------------	--------------------------

Returns

Returns the filtered second encoder rotational speed value

5.10.2.16 filter_vel_3()

```
int32 filter_vel_3 (  
    int32 value )
```

Filter on third encoder rotational speed (if present). The weighted average between the old value and the new one is executed.

Parameters

<i>value</i>	New value of the filter.
--------------	--------------------------

Returns

Returns the filtered third encoder rotational speed value

5.10.2.17 filter_voltage()

```
int32 filter_voltage (
    int32 value )
```

Filter on voltage readings. The weighted average between the old value and the new one is executed.

Parameters

<i>value</i>	New value of the filter.
--------------	--------------------------

Returns

Returns the filtered voltage value

5.10.2.18 my_mod()

```
uint32 my_mod (
    int32 val,
    int32 divisor )
```

This function computes the module function, returning positive values regardless of wheter the value passed is negative

Parameters

<i>val</i>	The value of which the module needs to be calculated
<i>divisor</i>	The divisor according to which the module is calculated

5.10.2.19 my_round()

```
int my_round (
    const double x )
```

This functions approximates the value passed to the nearest integer

Parameters

<i>x</i>	The floating point value that needs to be rounded
----------	---

Index

- ANTI_WINDUP
 - globals.h, 46
- acc
 - st_meas, 10
- activ
 - st_mem, 12
- activate_pwm_rescaling
 - st_mem, 13
- analog_read_end
 - interruptions.c, 54
 - interruptions.h, 58
- baud_rate
 - st_mem, 13
- buffer
 - st_data, 8
- c_mem
 - globals.c, 40
 - globals.h, 50
- CALIBRATION_DIV
 - globals.h, 46
- CURR_INTEGRAL_SAT_LIMIT
 - globals.h, 47
- CURRENT_HYSTERESIS
 - globals.h, 47
- CY_ISR_PROTO
 - interruptions.h, 58
- calc_turns_fcn
 - utils.c, 63
 - utils.h, 71
- calibration
 - utils.c, 63
 - utils.h, 71
- check_enc_data
 - utils.c, 63
 - utils.h, 71
- check_rest_position
 - utils.c, 63
 - utils.h, 72
- cmd_activate
 - command_processing.c, 21
 - command_processing.h, 29
- cmd_get_activate
 - command_processing.c, 21
 - command_processing.h, 29
- cmd_get_currents
 - command_processing.c, 21
 - command_processing.h, 29
- cmd_get_emg
 - command_processing.c, 21
 - command_processing.h, 29
- cmd_get_inputs
 - command_processing.c, 21
 - command_processing.h, 29
- cmd_get_measurements
 - command_processing.c, 21
 - command_processing.h, 29
- cmd_get_measurements_from_SH
 - command_processing.c, 21
 - command_processing.h, 29
- cmd_ping
 - command_processing.c, 21
 - command_processing.h, 29
- cmd_set_baudrate
 - command_processing.c, 22
 - command_processing.h, 30
- cmd_set_inputs
 - command_processing.c, 22
 - command_processing.h, 30
- cmd_set_watchdog
 - command_processing.c, 22
 - command_processing.h, 30
- cmd_store_params
 - command_processing.c, 22
 - command_processing.h, 30
- commProcess
 - command_processing.c, 22
 - command_processing.h, 30
- commReadMeasFromAnother
 - command_processing.c, 22
 - command_processing.h, 30
- commWrite
 - command_processing.c, 22
 - command_processing.h, 30
- commWrite_old_id
 - command_processing.c, 23
 - command_processing.h, 31
- commWriteAnother
 - command_processing.c, 23
 - command_processing.h, 31
- command_processing.c, 19
 - cmd_activate, 21
 - cmd_get_activate, 21
 - cmd_get_currents, 21
 - cmd_get_emg, 21
 - cmd_get_inputs, 21
 - cmd_get_measurements, 21
 - cmd_get_measurements_from_SH, 21

- cmd_ping, 21
- cmd_set_baudrate, 22
- cmd_set_inputs, 22
- cmd_set_watchdog, 22
- cmd_store_params, 22
- commProcess, 22
- commReadMeasFromAnother, 22
- commWrite, 22
- commWrite_old_id, 23
- commWriteAnother, 23
- get_param_list, 23
- infoGet, 24
- infoPrepare, 24
- infoSend, 24
- LCRChecksum, 24
- memInit, 25
- memRecall, 25
- memRestore, 25
- memStore, 25
- sendAcknowledgment, 26
- setZeros, 26
- command_processing.h, 26
 - cmd_activate, 29
 - cmd_get_activate, 29
 - cmd_get_currents, 29
 - cmd_get_emg, 29
 - cmd_get_inputs, 29
 - cmd_get_measurements, 29
 - cmd_get_measurements_from_SH, 29
 - cmd_ping, 29
 - cmd_set_baudrate, 30
 - cmd_set_inputs, 30
 - cmd_set_watchdog, 30
 - cmd_store_params, 30
 - commProcess, 30
 - commReadMeasFromAnother, 30
 - commWrite, 30
 - commWrite_old_id, 31
 - commWriteAnother, 31
 - get_param_list, 31
 - infoGet, 32
 - infoPrepare, 32
 - infoSend, 32
 - LCRChecksum, 32
 - memInit, 33
 - memRecall, 33
 - memRestore, 33
 - memStore, 33
 - sendAcknowledgment, 34
 - setZeros, 34
- command_slave
 - interruptions.c, 54
 - interruptions.h, 58
- commands.h, 34
 - qbmove_command, 36
 - qbmove_control_mode, 38
 - qbmove_input_mode, 38
 - qbmove_parameter, 38
- control_mode
 - st_mem, 13
- curr
 - st_meas, 11
 - st_ref, 18
- curr_estim
 - utils.c, 64
 - utils.h, 72
- curr_lookup
 - st_mem, 13
- current_limit
 - st_mem, 13
- DEFAULT_CURRENT_LIMIT
 - globals.h, 47
- DEFAULT_EEPROM_DISPLACEMENT
 - globals.h, 47
- DIV_INIT_VALUE
 - globals.h, 47
- deactivate_slave
 - interruptions.c, 54
 - interruptions.h, 58
- dev_pwm_limit
 - globals.c, 41
 - globals.h, 50
- dev_tension
 - globals.c, 41
 - globals.h, 50
- dev_tension_f
 - globals.c, 41
 - globals.h, 50
- direction
 - st_calib, 7
- double_encoder_on_off
 - st_mem, 13
- EMG_SAMPLE_TO_DISCARD
 - globals.h, 47
- emg
 - st_meas, 11
- emg_calibration_flag
 - st_mem, 13
- emg_max_value
 - st_mem, 13
- emg_speed
 - st_mem, 14
- emg_status
 - globals.h, 49
- emg_threshold
 - st_mem, 14
- enabled
 - st_calib, 7
- encoder_reading
 - interruptions.c, 54
 - interruptions.h, 59
- filter_acc_1
 - utils.c, 64
 - utils.h, 72

- filter_acc_2
 - utils.c, 64
 - utils.h, 73
- filter_acc_3
 - utils.c, 65
 - utils.h, 73
- filter_ch1
 - utils.c, 65
 - utils.h, 73
- filter_ch2
 - utils.c, 65
 - utils.h, 74
- filter_curr_diff
 - utils.c, 66
 - utils.h, 74
- filter_i1
 - utils.c, 66
 - utils.h, 75
- filter_v
 - utils.c, 67
 - utils.h, 75
- filter_vel_1
 - utils.c, 67
 - utils.h, 75
- filter_vel_2
 - utils.c, 67
 - utils.h, 76
- filter_vel_3
 - utils.c, 68
 - utils.h, 76
- filter_voltage
 - utils.c, 68
 - utils.h, 76
- flag
 - st_mem, 14
- forced_open
 - globals.c, 41
 - globals.h, 50
- function_scheduler
 - interruptions.c, 55
 - interruptions.h, 59
- g_measOld
 - globals.c, 41
 - globals.h, 50
- g_refOld
 - globals.c, 41
 - globals.h, 51
- g_rx
 - globals.c, 41
 - globals.h, 51
- get_param_list
 - command_processing.c, 23
 - command_processing.h, 31
- globals.c, 39
 - c_mem, 40
 - dev_pwm_limit, 41
 - dev_tension, 41
 - dev_tension_f, 41
 - forced_open, 41
 - g_measOld, 41
 - g_refOld, 41
 - g_rx, 41
 - interrupt_flag, 41
 - master_mode, 42
 - pow_tension, 42
 - pwm_sign, 42
 - receive_meas_from_hand, 42
 - reset_last_value_flag, 42
 - rest_enabled, 42
 - rest_pos_curr_ref, 42
 - tau_feedback, 42
 - tension_valid, 43
 - timer_value, 43
 - timer_value0, 43
 - watchdog_flag, 43
- globals.h, 43
 - ANTI_WINDUP, 46
 - c_mem, 50
 - CALIBRATION_DIV, 46
 - CURR_INTEGRAL_SAT_LIMIT, 47
 - CURRENT_HYSTERESIS, 47
 - DEFAULT_CURRENT_LIMIT, 47
 - DEFAULT_EEPROM_DISPLACEMENT, 47
 - DIV_INIT_VALUE, 47
 - dev_pwm_limit, 50
 - dev_tension, 50
 - dev_tension_f, 50
 - EMG_SAMPLE_TO_DISCARD, 47
 - emg_status, 49
 - forced_open, 50
 - g_measOld, 50
 - g_refOld, 51
 - g_rx, 51
 - interrupt_flag, 51
 - LOOKUP_DIM, 47
 - MAX_WATCHDOG_TIMER, 47
 - master_mode, 51
 - NUM_OF_ANALOG_INPUTS, 48
 - NUM_OF_EMGS, 48
 - NUM_OF_MOTORS, 48
 - NUM_OF_PARAMS, 48
 - NUM_OF_SENSORS, 48
 - POS_INTEGRAL_SAT_LIMIT, 48
 - PWM_MAX_VALUE, 48
 - pow_tension, 51
 - pwm_sign, 51
 - RECEIVE, 48
 - receive_meas_from_hand, 51
 - reset_last_value_flag, 51
 - rest_enabled, 52
 - rest_pos_curr_ref, 52
 - SAMPLES_FOR_EMG_MEAN, 49
 - SAMPLES_FOR_MEAN, 49
 - tau_feedback, 52
 - tension_valid, 52
 - timer_value, 52

- timer_value0, 52
- UNLOAD, 49
- WAIT_ID, 49
- WAIT_LENGTH, 49
- WAIT_START, 49
- watchdog_flag, 52
- hand_meas
 - st_meas, 11
- id
 - st_mem, 14
- ind
 - st_data, 8
- infoGet
 - command_processing.c, 24
 - command_processing.h, 32
- infoPrepare
 - command_processing.c, 24
 - command_processing.h, 32
- infoSend
 - command_processing.c, 24
 - command_processing.h, 32
- input_mode
 - st_mem, 14
- interrupt_flag
 - globals.c, 41
 - globals.h, 51
- interrupt_manager
 - interruptions.c, 55
 - interruptions.h, 59
- interruptions.c, 53
 - analog_read_end, 54
 - command_slave, 54
 - deactivate_slave, 54
 - encoder_reading, 54
 - function_scheduler, 55
 - interrupt_manager, 55
 - motor_control, 55
 - overcurrent_control, 55
 - pwm_limit_search, 55
 - pwm_preload_values, 55
- interruptions.h, 56
 - analog_read_end, 58
 - CY_ISR_PROTO, 58
 - command_slave, 58
 - deactivate_slave, 58
 - encoder_reading, 59
 - function_scheduler, 59
 - interrupt_manager, 59
 - motor_control, 59
 - overcurrent_control, 59
 - pwm_limit_search, 59
- k_d
 - st_mem, 14
- k_d_c
 - st_mem, 14
- k_d_c_dl
 - st_mem, 14
- k_d_dl
 - st_mem, 15
- k_i
 - st_mem, 15
- k_i_c
 - st_mem, 15
- k_i_c_dl
 - st_mem, 15
- k_i_dl
 - st_mem, 15
- k_p
 - st_mem, 15
- k_p_c
 - st_mem, 15
- k_p_c_dl
 - st_mem, 15
- k_p_dl
 - st_mem, 16
- LCRChecksum
 - command_processing.c, 24
 - command_processing.h, 32
- LOOKUP_DIM
 - globals.h, 47
- length
 - st_data, 9
- m_mult
 - st_mem, 16
- m_off
 - st_mem, 16
- MAX_WATCHDOG_TIMER
 - globals.h, 47
- main.c, 60
- master_mode
 - globals.c, 42
 - globals.h, 51
- max_step_neg
 - st_mem, 16
- max_step_pos
 - st_mem, 16
- memInit
 - command_processing.c, 25
 - command_processing.h, 33
- memRecall
 - command_processing.c, 25
 - command_processing.h, 33
- memRestore
 - command_processing.c, 25
 - command_processing.h, 33
- memStore
 - command_processing.c, 25
 - command_processing.h, 33
- motor_control
 - interruptions.c, 55
 - interruptions.h, 59
- motor_handle_ratio
 - st_mem, 16

- my_mod
 - utils.c, 68
 - utils.h, 77
- my_round
 - utils.c, 69
 - utils.h, 77
- NUM_OF_ANALOG_INPUTS
 - globals.h, 48
- NUM_OF_EMGS
 - globals.h, 48
- NUM_OF_MOTORS
 - globals.h, 48
- NUM_OF_PARAMS
 - globals.h, 48
- NUM_OF_SENSORS
 - globals.h, 48
- onoff
 - st_ref, 18
- overcurrent_control
 - interruptions.c, 55
 - interruptions.h, 59
- POS_INTEGRAL_SAT_LIMIT
 - globals.h, 48
- PWM_MAX_VALUE
 - globals.h, 48
- pos
 - st_meas, 11
 - st_ref, 18
- pos_lim_flag
 - st_mem, 16
- pos_lim_inf
 - st_mem, 16
- pos_lim_sup
 - st_mem, 17
- pow_tension
 - globals.c, 42
 - globals.h, 51
- pwm
 - st_ref, 18
- pwm_limit
 - st_dev, 9
- pwm_limit_search
 - interruptions.c, 55
 - interruptions.h, 59
- pwm_preload_values
 - interruptions.c, 55
- pwm_sign
 - globals.c, 42
 - globals.h, 51
- qbmove_command
 - commands.h, 36
- qbmove_control_mode
 - commands.h, 38
- qbmove_input_mode
 - commands.h, 38
- qbmove_parameter
 - commands.h, 38
- RECEIVE
 - globals.h, 48
- ready
 - st_data, 9
- receive_meas_from_hand
 - globals.c, 42
 - globals.h, 51
- repetitions
 - st_calib, 7
- res
 - st_mem, 17
- reset_last_value_flag
 - globals.c, 42
 - globals.h, 51
- rest_delay
 - st_mem, 17
- rest_enabled
 - globals.c, 42
 - globals.h, 52
- rest_pos
 - st_mem, 17
- rest_pos_curr_ref
 - globals.c, 42
 - globals.h, 52
- rest_ratio
 - st_mem, 17
- rest_vel
 - st_mem, 17
- rot
 - st_meas, 11
- SAMPLES_FOR_EMG_MEAN
 - globals.h, 49
- SAMPLES_FOR_MEAN
 - globals.h, 49
- sendAcknowledgment
 - command_processing.c, 26
 - command_processing.h, 34
- setZeros
 - command_processing.c, 26
 - command_processing.h, 34
- speed
 - st_calib, 8
- st_calib, 7
 - direction, 7
 - enabled, 7
 - repetitions, 7
 - speed, 8
- st_data, 8
 - buffer, 8
 - ind, 8
 - length, 9
 - ready, 9
- st_dev, 9
 - pwm_limit, 9
 - tension, 10

- tension_conv_factor, 10
- tension_valid, 10
- st_meas, 10
 - acc, 10
 - curr, 11
 - emg, 11
 - hand_meas, 11
 - pos, 11
 - rot, 11
 - vel, 11
- st_mem, 12
 - activ, 12
 - activate_pwm_rescaling, 13
 - baud_rate, 13
 - control_mode, 13
 - curr_lookup, 13
 - current_limit, 13
 - double_encoder_on_off, 13
 - emg_calibration_flag, 13
 - emg_max_value, 13
 - emg_speed, 14
 - emg_threshold, 14
 - flag, 14
 - id, 14
 - input_mode, 14
 - k_d, 14
 - k_d_c, 14
 - k_d_c_dl, 14
 - k_d_dl, 15
 - k_i, 15
 - k_i_c, 15
 - k_i_c_dl, 15
 - k_i_dl, 15
 - k_p, 15
 - k_p_c, 15
 - k_p_c_dl, 15
 - k_p_dl, 16
 - m_mult, 16
 - m_off, 16
 - max_step_neg, 16
 - max_step_pos, 16
 - motor_handle_ratio, 16
 - pos_lim_flag, 16
 - pos_lim_inf, 16
 - pos_lim_sup, 17
 - res, 17
 - rest_delay, 17
 - rest_pos, 17
 - rest_ratio, 17
 - rest_vel, 17
 - watchdog_period, 17
- st_ref, 18
 - curr, 18
 - onoff, 18
 - pos, 18
 - pwm, 18
- tau_feedback
 - globals.c, 42
 - globals.h, 52
- tension
 - st_dev, 10
- tension_conv_factor
 - st_dev, 10
- tension_valid
 - globals.c, 43
 - globals.h, 52
 - st_dev, 10
- timer_value
 - globals.c, 43
 - globals.h, 52
- timer_value0
 - globals.c, 43
 - globals.h, 52
- UNLOAD
 - globals.h, 49
- utils.c, 61
 - calc_turns_fcn, 63
 - calibration, 63
 - check_enc_data, 63
 - check_rest_position, 63
 - curr_estim, 64
 - filter_acc_1, 64
 - filter_acc_2, 64
 - filter_acc_3, 65
 - filter_ch1, 65
 - filter_ch2, 65
 - filter_curr_diff, 66
 - filter_i1, 66
 - filter_v, 67
 - filter_vel_1, 67
 - filter_vel_2, 67
 - filter_vel_3, 68
 - filter_voltage, 68
 - my_mod, 68
 - my_round, 69
- utils.h, 69
 - calc_turns_fcn, 71
 - calibration, 71
 - check_enc_data, 71
 - check_rest_position, 72
 - curr_estim, 72
 - filter_acc_1, 72
 - filter_acc_2, 73
 - filter_acc_3, 73
 - filter_ch1, 73
 - filter_ch2, 74
 - filter_curr_diff, 74
 - filter_i1, 75
 - filter_v, 75
 - filter_vel_1, 75
 - filter_vel_2, 76
 - filter_vel_3, 76
 - filter_voltage, 76
 - my_mod, 77
 - my_round, 77

vel
 st_meas, 11

WAIT_ID
 globals.h, 49

WAIT_LENGTH
 globals.h, 49

WAIT_START
 globals.h, 49

watchdog_flag
 globals.c, 43
 globals.h, 52

watchdog_period
 st_mem, 17